


### Monolayer™/Grazing Angle Accessory for FT-IR Spectroscopy

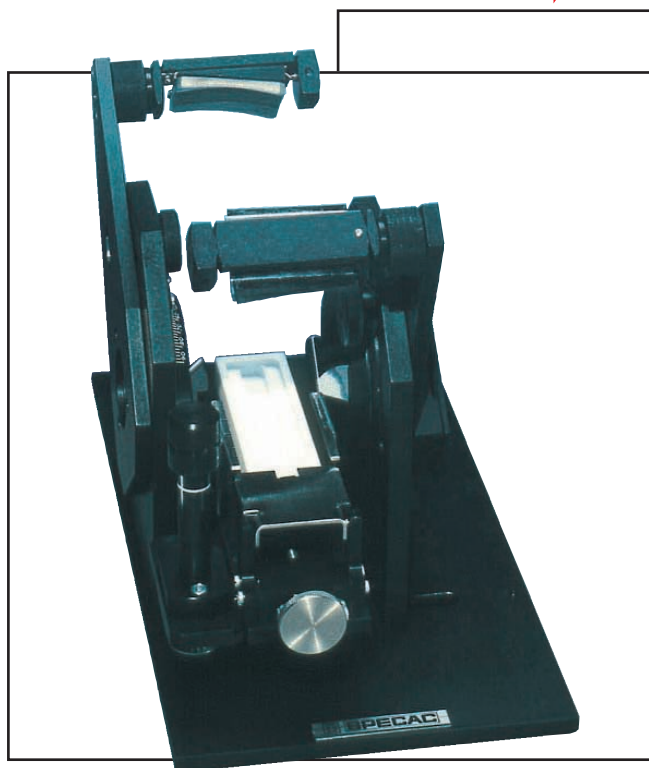

**A grazing incidence specular reflectance accessory operating at ambient temperatures.**

#### Features

- Simple alignment
- Easy accessory interchangeability
- High throughput
- Wide incidence angle 8°- 85°
- Interchangeable aperture stops
- Polarizer mount
- Film expansion/compression facility

The Monolayer™/ Grazing Angle Accessory allows the study of insoluble monolayer films at air-water interface in situ using an MCT or InSb detector, and also the analysis of reflective surfaces at grazing incidence angle (8°-85°) using a standard DTGS detector.

Alignment is very easy, as it involves the adjustments of only two mirrors. Conversion from one form to another is achieved by a simple change of sample holders.

In the monolayer mode, a variety of monolayer systems (eg. Surfactants, proteins, steroids, polymers and phospholipids in biological membranes) can be studied in their native environments. The ability to stretch and compress the films is incorporated, and useful information such as phase transition changes, solvation of head group, orientation of hydrophobic chain, conformation, configuration and polarization studies may be made.

Operating in the grazing incidence mode at high angles of incidence, up to 85°, allows analysis of thin layers of Langmuir-Blodgett films adsorbed on reflective surfaces. Aperture stops are provided to redefine the beam profile to be at prime position on samples at or near grazing angle to enhance the performance. In this operation mode involving air-metal interface, optimization of spectral absorption of coatings on reflective backgrounds and performance data of semi-conductor materials and films can be evaluated.

#### Ordering Information

19650 Monolayer™/Grazing Angle Accessory with aperture stops (please specify FT-IR make and model).

*\*\* A wide range of polarizers is offered (details on request).*

THE SOLUTION IS SPECAC WHAT'S THE PROBLEM?

smiths  
A part of Smiths Group plc

# Specac

www.specac.com



**SPECAC LIMITED,**  
River House, 97 Cray Ave, St. Mary Cray,  
Orpington, Kent BR5 4HE U.K.  
Tel : +44 (0)1689 873134  
Fax : +44 (0)1689 878527  
Registered No. 1008689 England

**SPECAC INC.,**  
410 Creekstone Ridge,  
Woodstock, GA 30188 U.S.A.  
Toll Free +1 800 447 2558  
www.specac.com  
Specac is part of Smiths Group plc