

**SAFETYCAPS ARE AVAILABLE WITH UP TO 6 TUBE CONNECTORS AND 6 SHUT-OFFS.**

## Removal systems for GL45 threaded bottles SafetyCap with shut-off

During repairs on HPLC pumps, the shut-off on the SafetyCaps can be closed to keep capillaries and inlet filters in the solvent. The frit does not come into contact with the ambient air. The pump can be quickly and easily flushed during purging after repairs. The SafetyCaps with shut-off also prevent air pockets in the tube when the supply bottle is being changed, and guarantee the ongoing analysis continues without a problem.



- different thread diameters
- 1 to 6 connectors
- also without shut-off → Page 10
- fire-resistant → Page 36
- optional: SafeLock stopper → Pages 147 + 26
- as set with containers → Page 18

Part No.	Description	Connectors
107119	SafetyCap I GL45	1 x connector Ø 3.2 mm OD with shut-off
107919	SafetyCap II GL45	2 x connector Ø 3.2 mm OD each with shut-off
107920	SafetyCap III GL45	3 x connector Ø 3.2 mm OD each with shut-off
107419	SafetyCap IV GL45	4 x connector Ø 3.2 mm OD each with shut-off <b>NEW</b>
107519	SafetyCap VI GL45	6 x connector Ø 3.2 mm OD each with shut-off <b>NEW</b>

All SafetyCaps are delivered with an air valve and corresponding fittings. Plugs and ferrules for capillaries with 1.6 or 2.3 mm outer diameter and other accessories can be found starting on → page 133.

### Accessories

Part No.	Description
117010	Air valve with air filter → Page 134

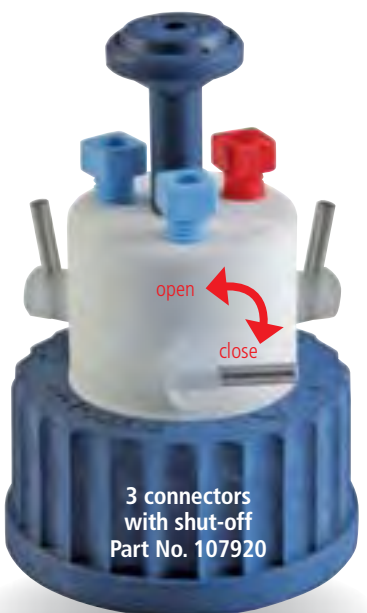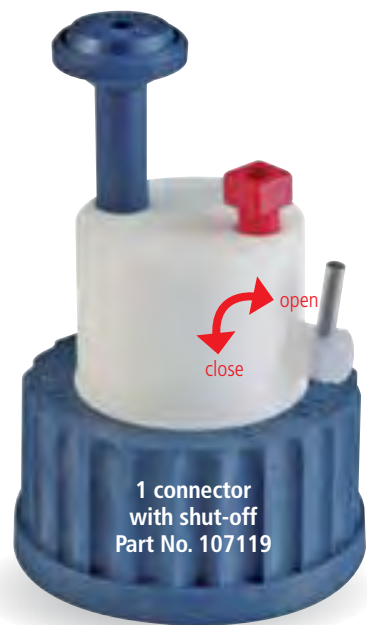

The shut-off allows a connector to be closed quickly and easily. This lets tubes and containers be changed securely and easily.

## Removal systems for GL45 threaded bottles SafetyCap combined

The SafetyCap Kombi and its all-around connector provide a stable baseline and reliable reproducibility for RI detection. For fast and simple solvent degassing with helium.

- different thread diameters
- 2 or 3 connectors
- fire-resistant → Page 36
- optional: SafeLock stopper → Pages 147 + 26
- as set with containers → Page 18

Part No.	Description	Connectors
107219	SafetyCap II (combined)	1x connector Ø 3.2 mm OD with shut-off and 1 x connector Ø 3.2 mm OD
107319	SafetyCap III (combined)	2 x connector Ø 3.2 mm OD with shut-off and 1 x connector Ø 3.2 mm OD

All SafetyCaps are delivered with an air valve and corresponding fittings. Plugs and ferrules for capillaries with 1.6 or 2.3 mm outer diameter and other accessories can be found starting on → page 133.

### Accessories

Part No.	Description
117010	Air valve with air filter → Page 134

The shut-off allows a connector to be closed quickly and easily. This lets tubes and containers be changed securely and easily.

