

Hydrocarbons C10 - C40

Column: **TRB-1HT**, P/N TR-610132

Size: 30m x 0.25mm x 0.10µm


Carrier gas: Hydrogen, 12psi (82.7kPa)

Sample: hydrocarbons mixture C10-C40(50µg/mL), C17,Pristane,C18,Phytane(30µg/mL)

Injection: 275°C, 1µL, split 1:25

Oven: 50°C to 350°C(5min) @ 10°C/min

Detector: FID,350°C

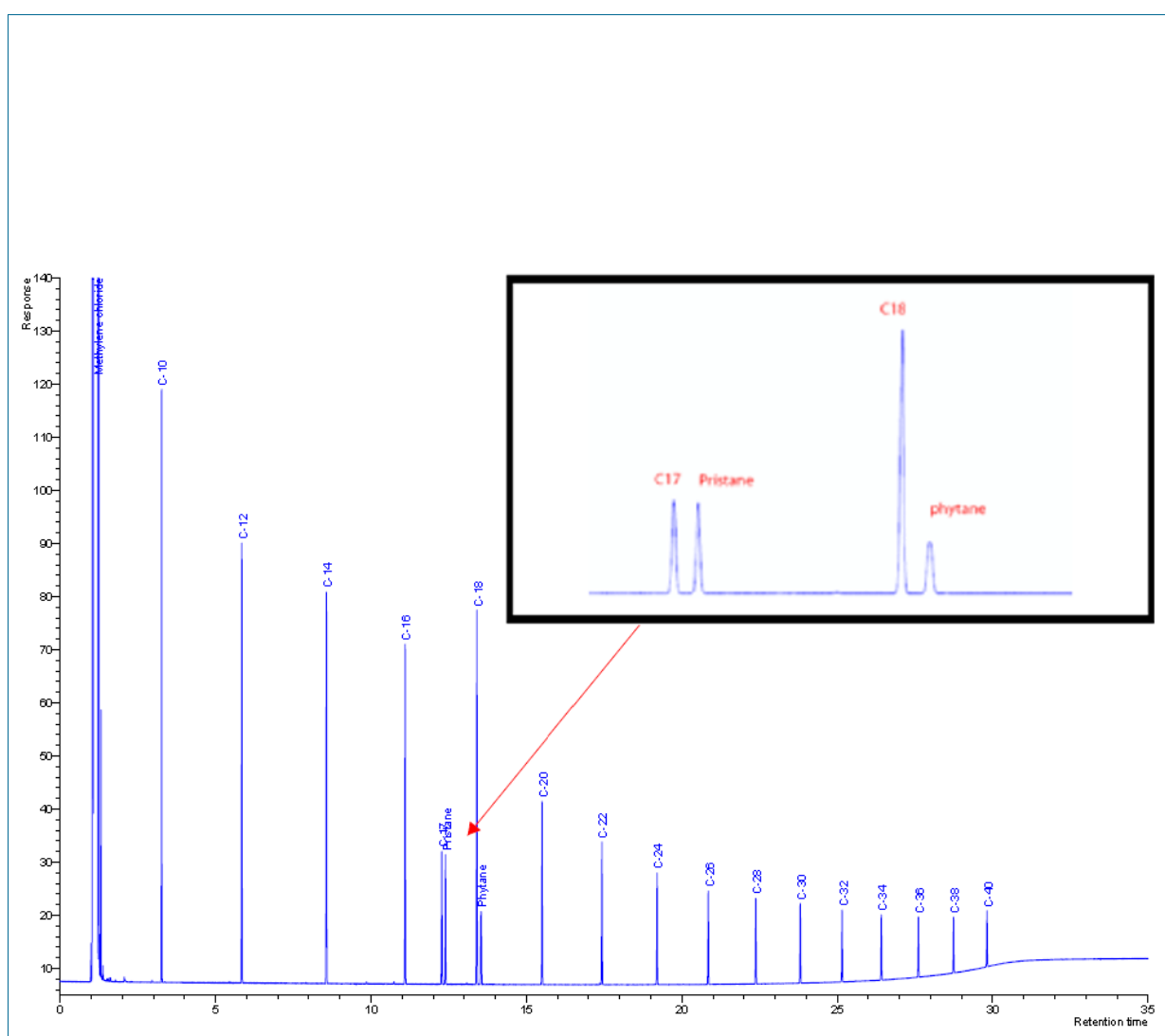

TKG 1283