

TRB-1ht

Column: Retention Gap (intermediate polarity) 5m x 0.53mm(TR-200055)

+**TRB-1ht** (TR-610113)


Dimensions: 15m x 0.32mm x 0.10µm

Injection: 0.3µL Poliwx 655 (0.1% in CS₂), on column (sec. cool 30s), 280°C

Carrier gas: H₂, constant flow 2mL/min

Oven temperature: 50°C(2 min) @ 15°C/min to 380°C(5 min)

Detector: FID, 390°C

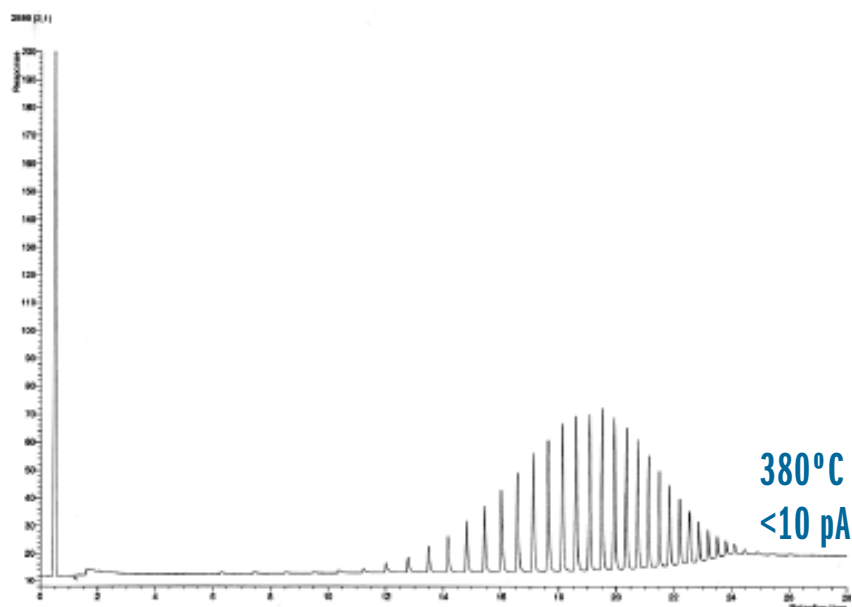

TKG 1103

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es