

HYDROCARBONS

Column: **TRB-1ht**, P/N TR-610133

Size: 30m x 0.32mm x 0.1µm


Injection: hydrocarbons standard 1250 ppb, splitless, 250°C

Carrier gas: He, constant flow 2 mL/min

Program Temperature: 50°C(1 min) @ 40°C/min to 320°C(5 min)

Detector: MS, ion source 200°C, Interface 280°C, scan 20-600

Chromatogram provided by Vanesa Riu de ILERSAP, Mollerussa (Lleida).

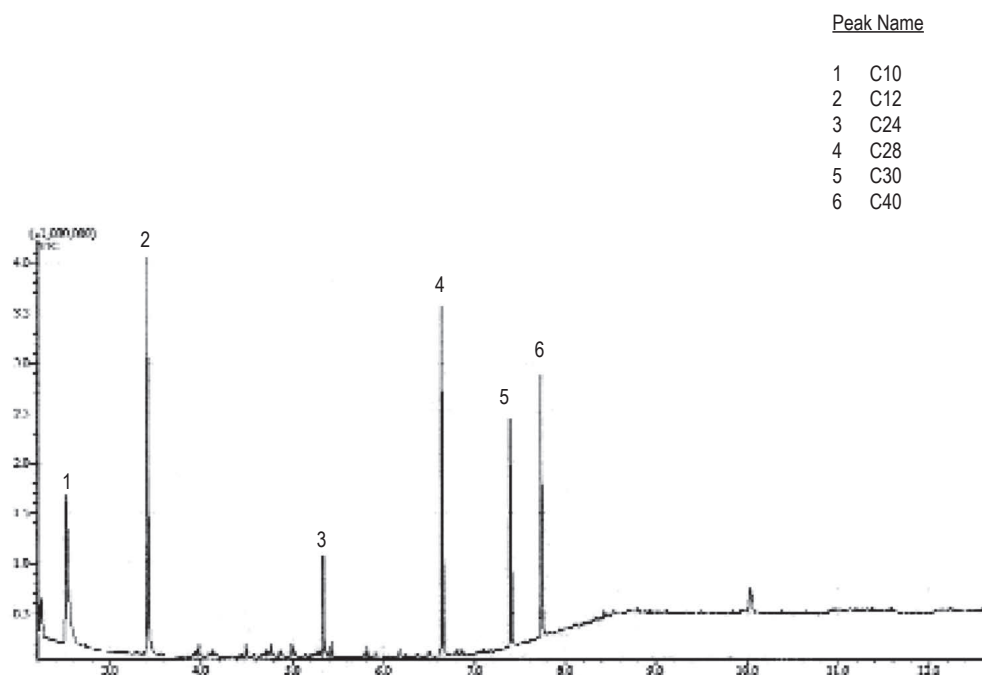

TKG 1221

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es