


TRB-2887

Column: **TBR-2887**, P/N TR-192645
 Dimensions: 10m x 0.53mm x 2.65µm
 Injection: 1µL standard RS-31222, split 1:10, 300°C
 Carrier gas: H₂, constant pressure 5 psi (34.4 KPa)
 Oven temperature: 40°C @ 15°C/min to 350°C(5 min)
 Detector: FID, 360°C. Baseline not compensated

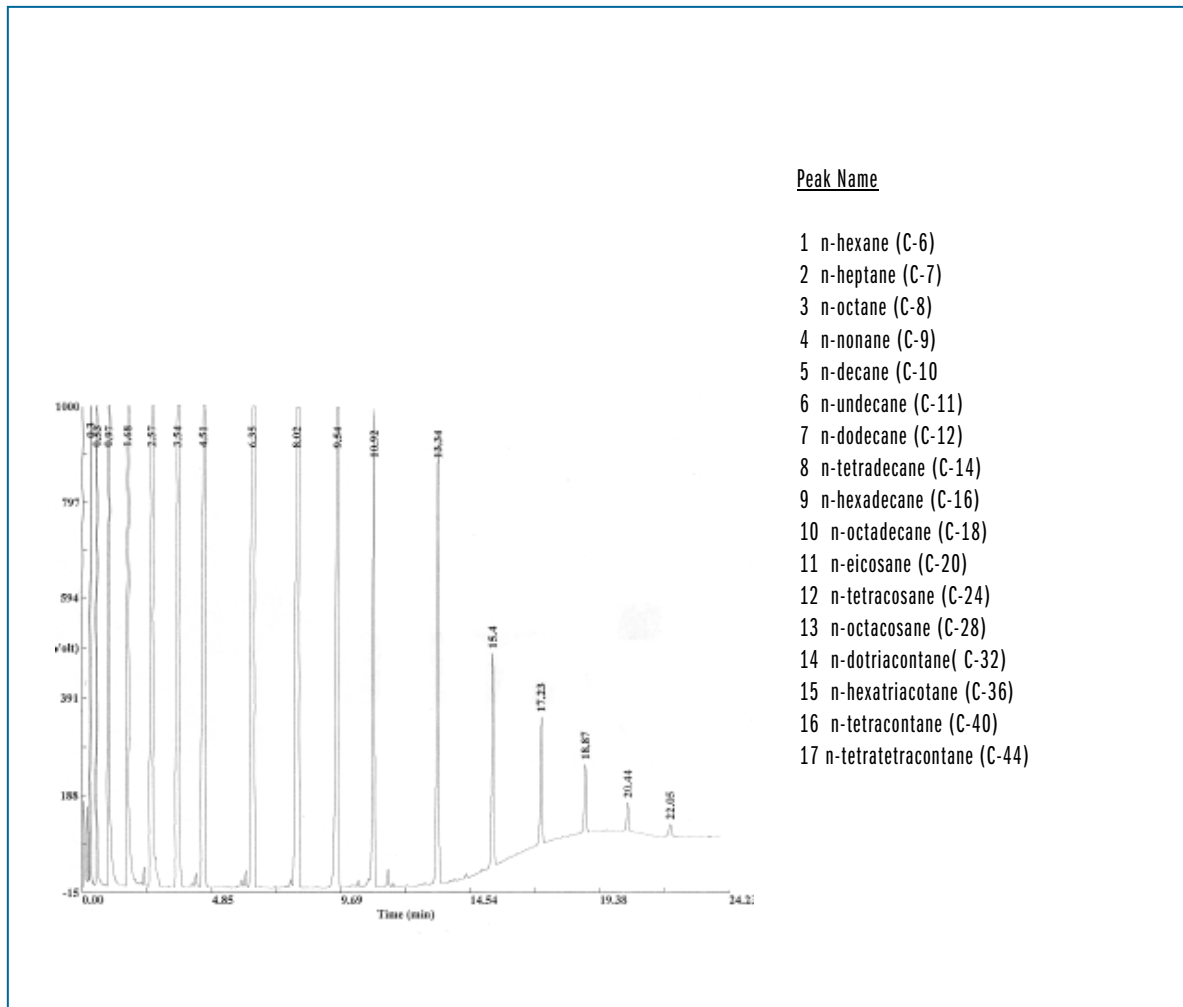

TKG 1108

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es