

STEROLS IN OLIVE OIL

Column: **TRB-STEROL**, P/N TR-180738

Size: 30m x 0.22mm x 0.12 μ m


Injection: 1 μ L olive oil in diethyl ether (extraction following norm), split 1:30, 250°C

Carrier Gas: H₂, constant flow 1.6 mL/min

Oven Temperature: 285°C

Detector: FID, 320°C

Chromatogram provided by Dr. Manuel León, Instituto de la Grasa (Sevilla).

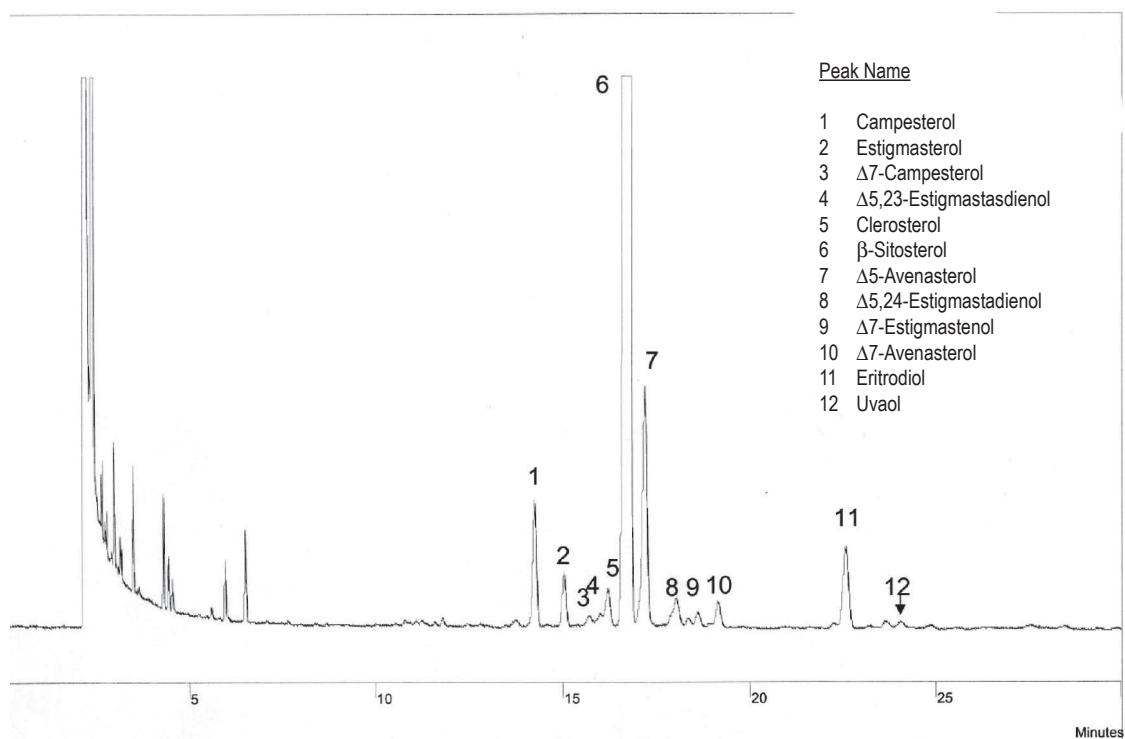

TKG 1197

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es