

Column: **MTI-5**, P/N TR-910232


Dimensions: 30m x 0.25mm x 0.25µm

Injection: 1µL Test MX5 (10 a 20 ng/comp), split 1:100, 280°C

Carrier Gas: H₂, 12 PSI (82.7 KPa)

Oven Temp: 100°C@ 6°C/min a 330° C (15min) @10°C/min a 360°C (15 min.)

Detector: FID, 360°C

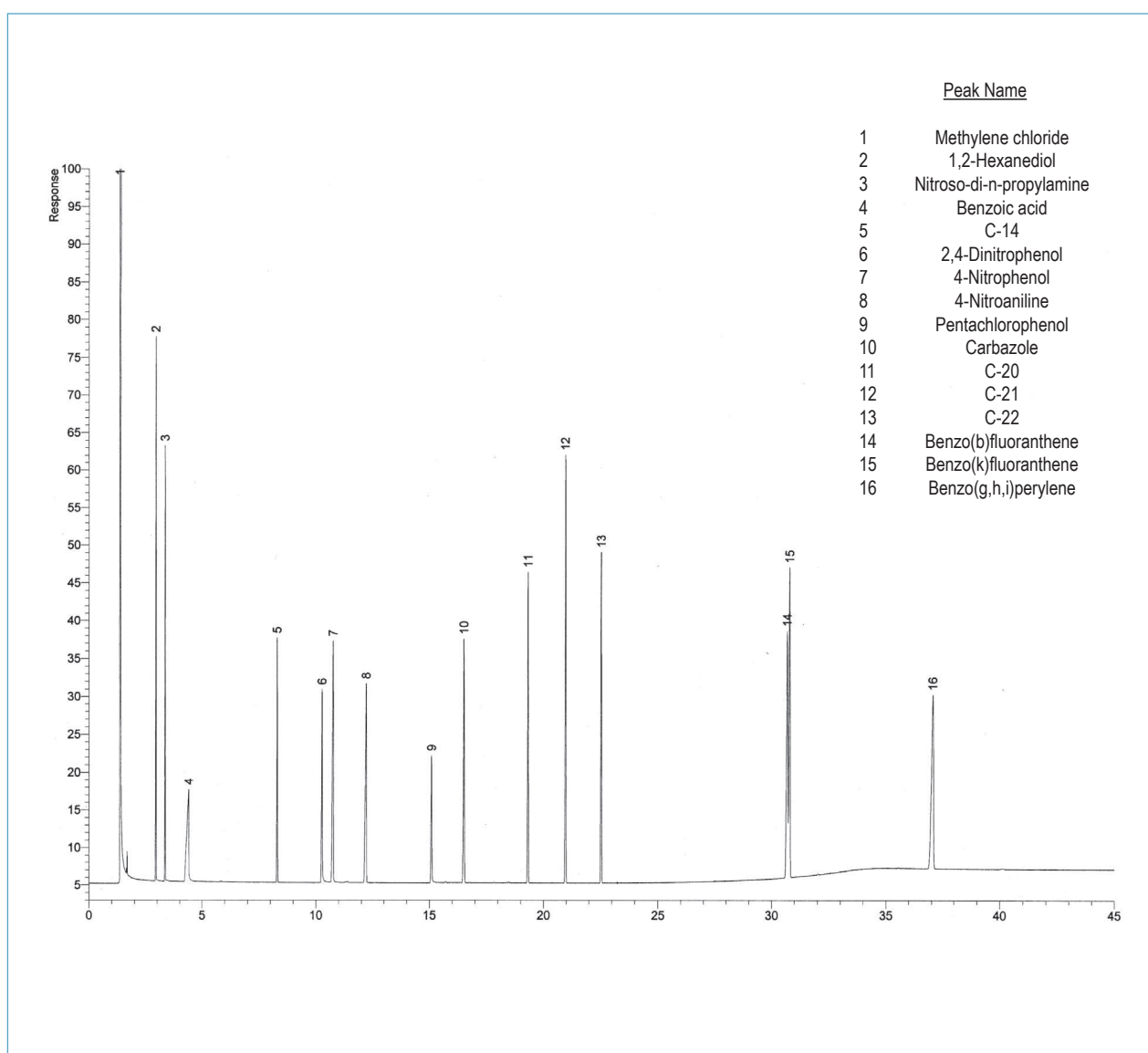

TKG 1207

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es