

TRB-1301

Column: **TRB-1301**, P/N TR-601032


Dimensions: 30m x 0.25mm x 1.0µm

Injection: 0.5µL standard SP-4-7301 (500 ng/mL), split 1:50, 260°C

Carrier gas: H₂, constant pressure 12 psi (82.7 KPa).

Oven temperature: 165°C

Detector: FID, 280°C

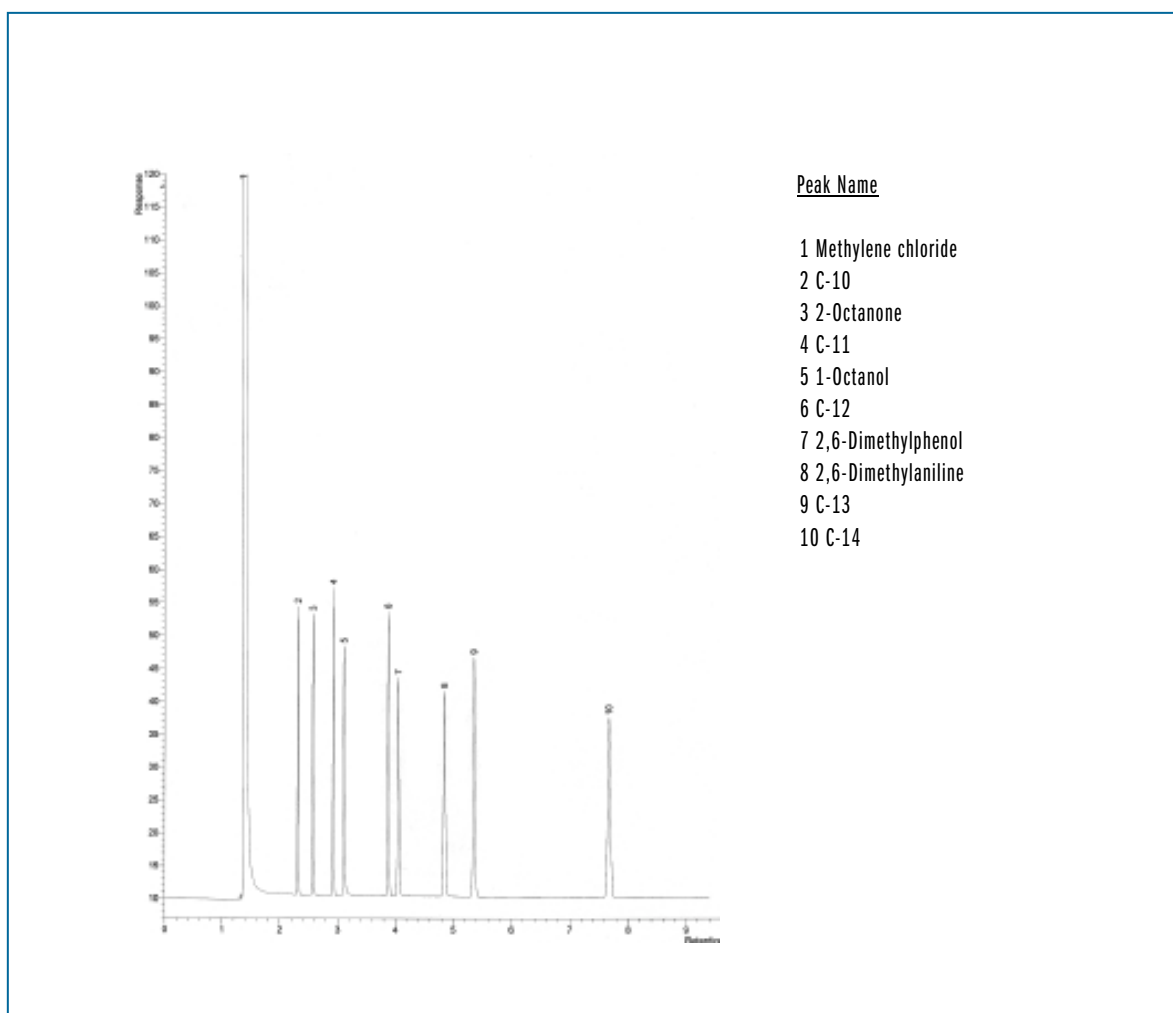

TKG 1118

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es