

TRB-G43

Column: **TRB-G43**, P/N TR-163035

Dimensions: 30m x 0.53mm x 3.0µm


Injection: split 1:2, 250°C, 5m x 0.53mm intermediate polarity column (TR-200055)

Sample: 0.5 mL headspace 80°C (2t static head space sampler) 28 Class 1 Mix
Class 2 Mix A, Mix B residual solvents at the regulatory limit concentration.

Carrier gas: He, constant pressure 4.8 psi (33.1 KPa), 35 cm/s (40°C)

Oven temperature: 40°C(20min) @ 10°C/min to 240°C(10min)

Detector: FID, 250°C

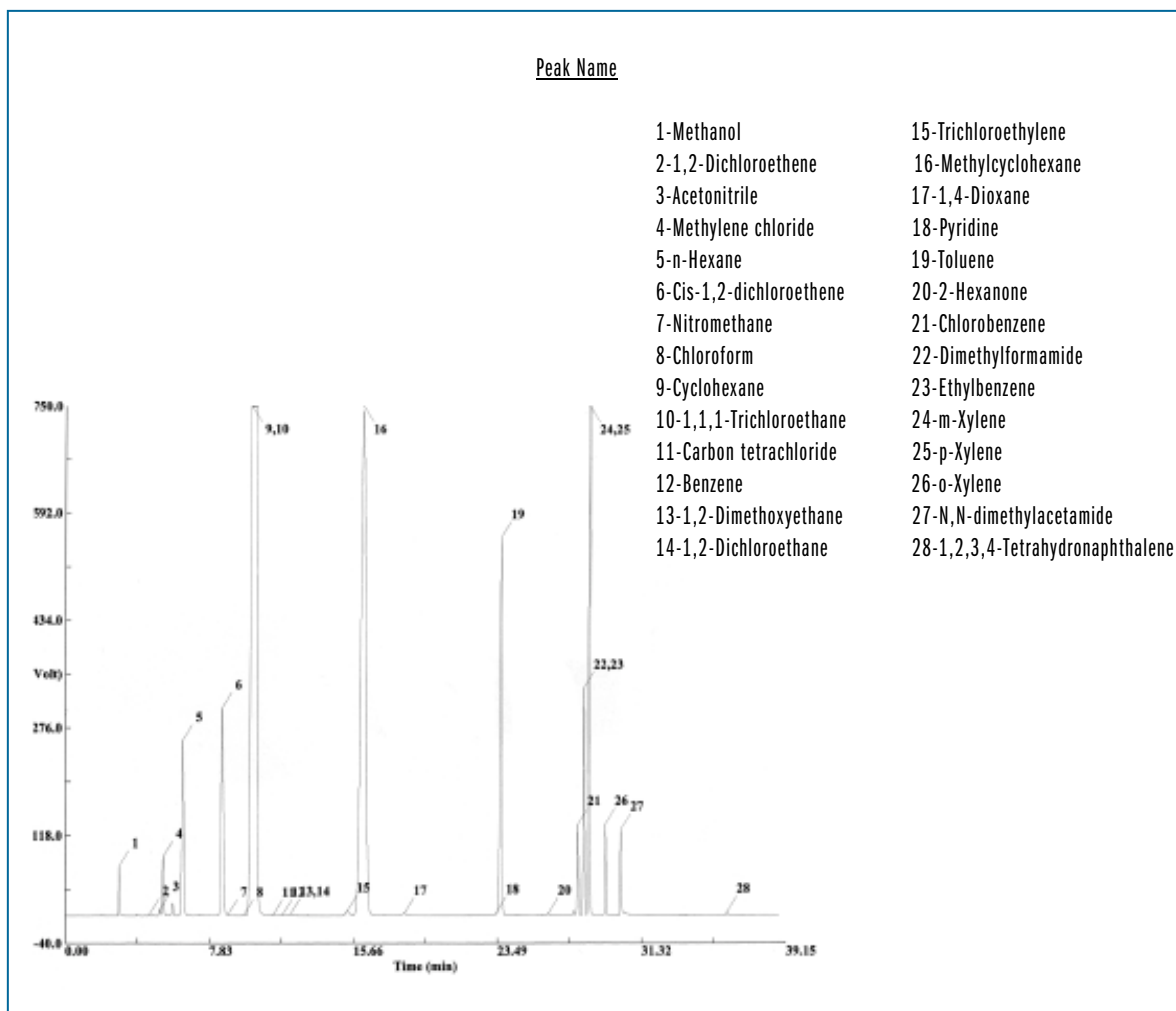

TKG 1121