


ALCOHOLS AND TERPENIC ACIDS (OLIVE OIL)

Column: **TRB-35**, P/N TR-351332
Dimensions: 30m x 0.25mm x 0.15µm
Injection: 1 µL extract of leaf of Olive Tree, split 1:20, 300°C
Carrier gas: H₂, constant pressure 12 psi (82.7 KPa).
Oven temperature: 275°C
Detector: FID, 300°C

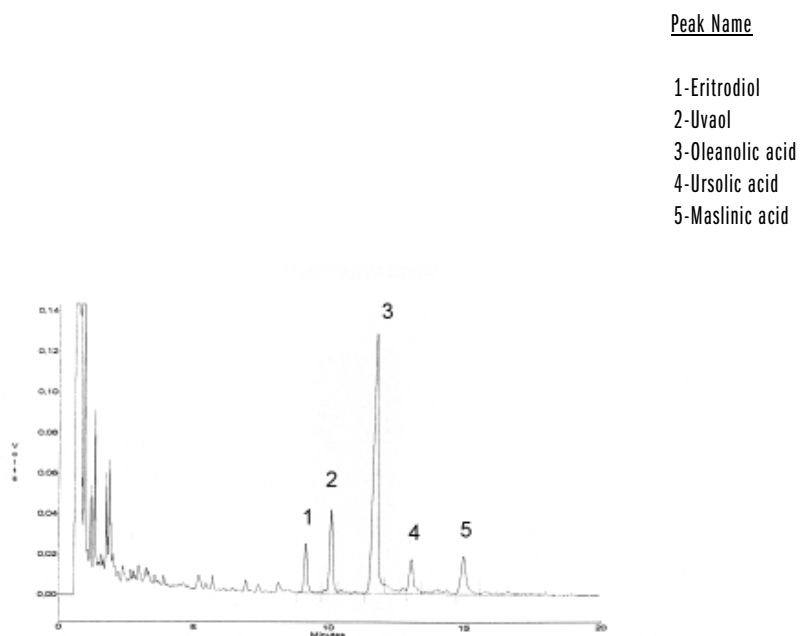

Chromatogram provided by Angeles Ginda from Instituto de la Grasa, CSIC.

TKG 1049