


FREE ACIDS IN WATER

Column: **TRB-FFAP**, P/N TR-151035
Dimensions: 30m x 0.53mm x 1.0µm
Injection: 1µL, split
Carrier gas: He, 4 psi (27.56 KPa)
Oven temperature: 120°C @ 4°C/min to 220°C
Detector: FID, 275°C

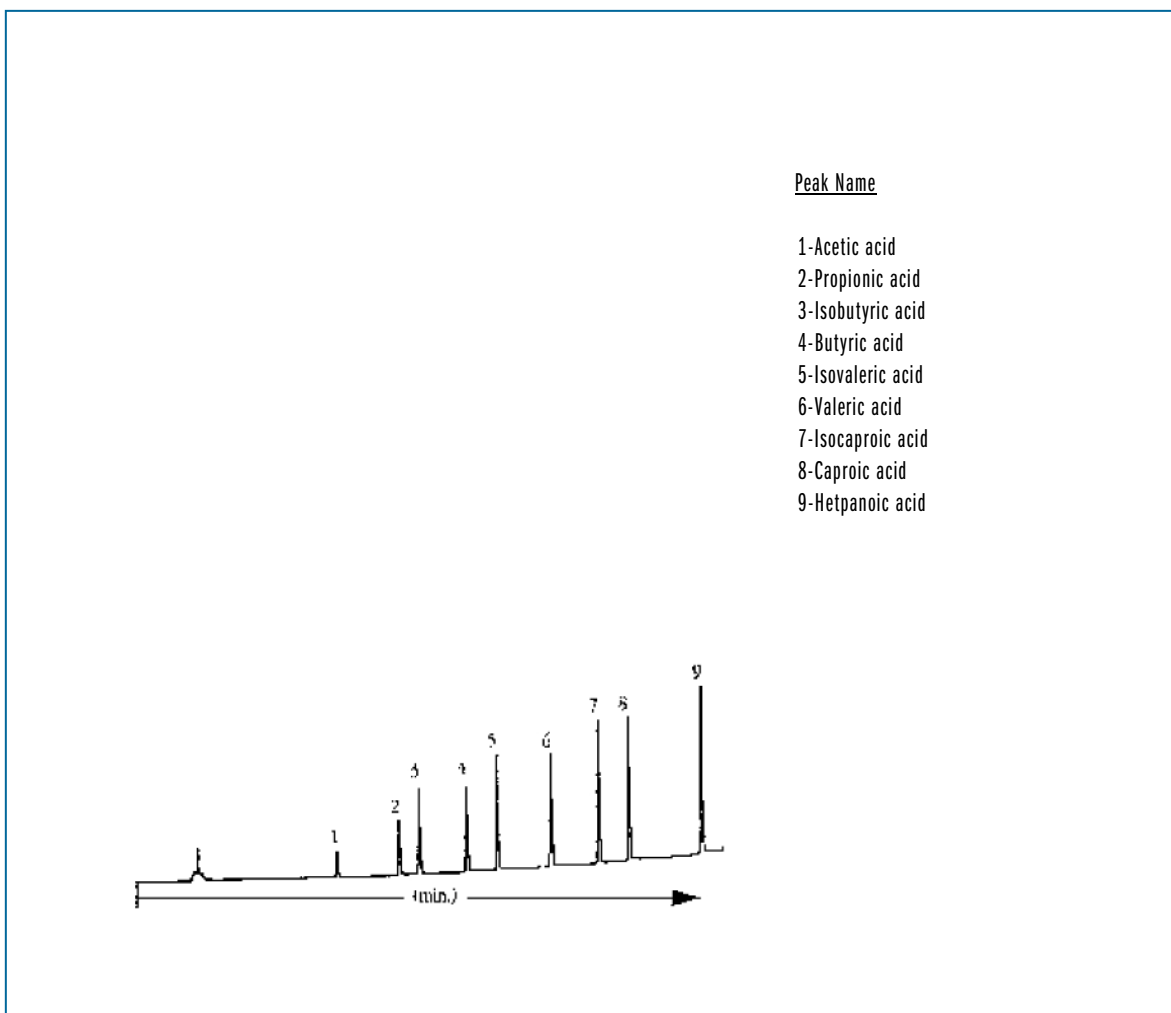

TKG 1148

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es