

ALCOHOLS IN BLOOD

Column: **MetaBLOOD 2**, P/N TR-862035

Size: 30m x 0.53mm x 2.0µm

Carrier gas: He, 5 psi

Oven Temperature: 45°C (15 min)


Injection: 1 mL Head Space 2t (vial 70°C), alcohols in blood (2-20 ppm), split 1:20, 225°C

Detector: FID, 300°C

*Chromatogram provided by Dra. Guadalupe Montoya y Dra. Isabel Bonaparte
from General Lab (Barcelona)*

Peak Name

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acetone
- 5 Isopropanol
- 6 n-Propanol
- 7 Methyl ethyl ketone (MEK)
- 8 Methylisobutylketone (MIBK)

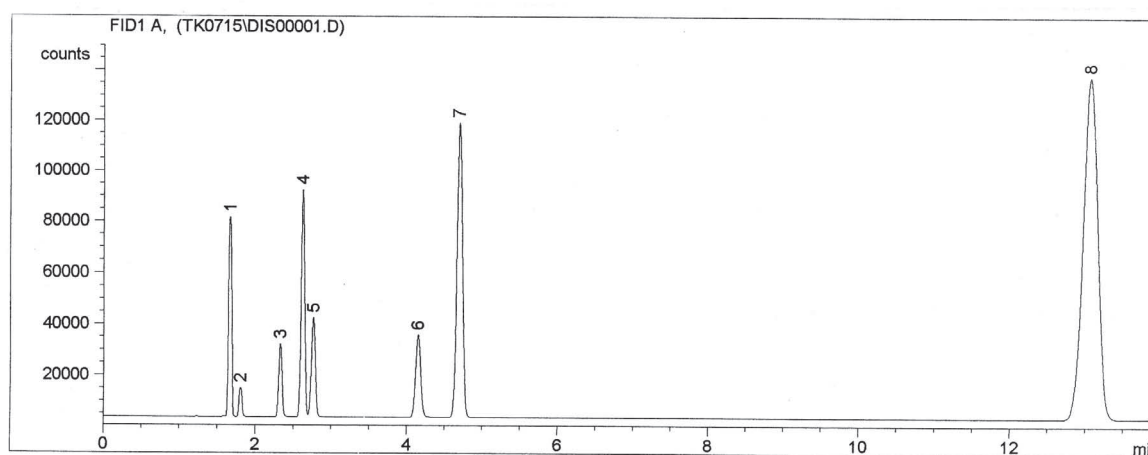

TKG 1209

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es