

SAPIENS-624MS: USP <467> Residual Solvents Class 2

Column: SAPIENS-624MS: P/N TR-591833


Dimensions: 30 m x 0,32 mm x 1,8 µm

Injection: 0.1 µl Residual Solvents Class 2 (0.98-19.4 mg/ml) in DMSO, 140°C, split: 1:50

Carrier gas: He, ct flow 2.2 ml/min

Oven temperature: 40°C(20 min) to 240°C (15 min) @ 10°C/min

Detector: FID @ 250°C

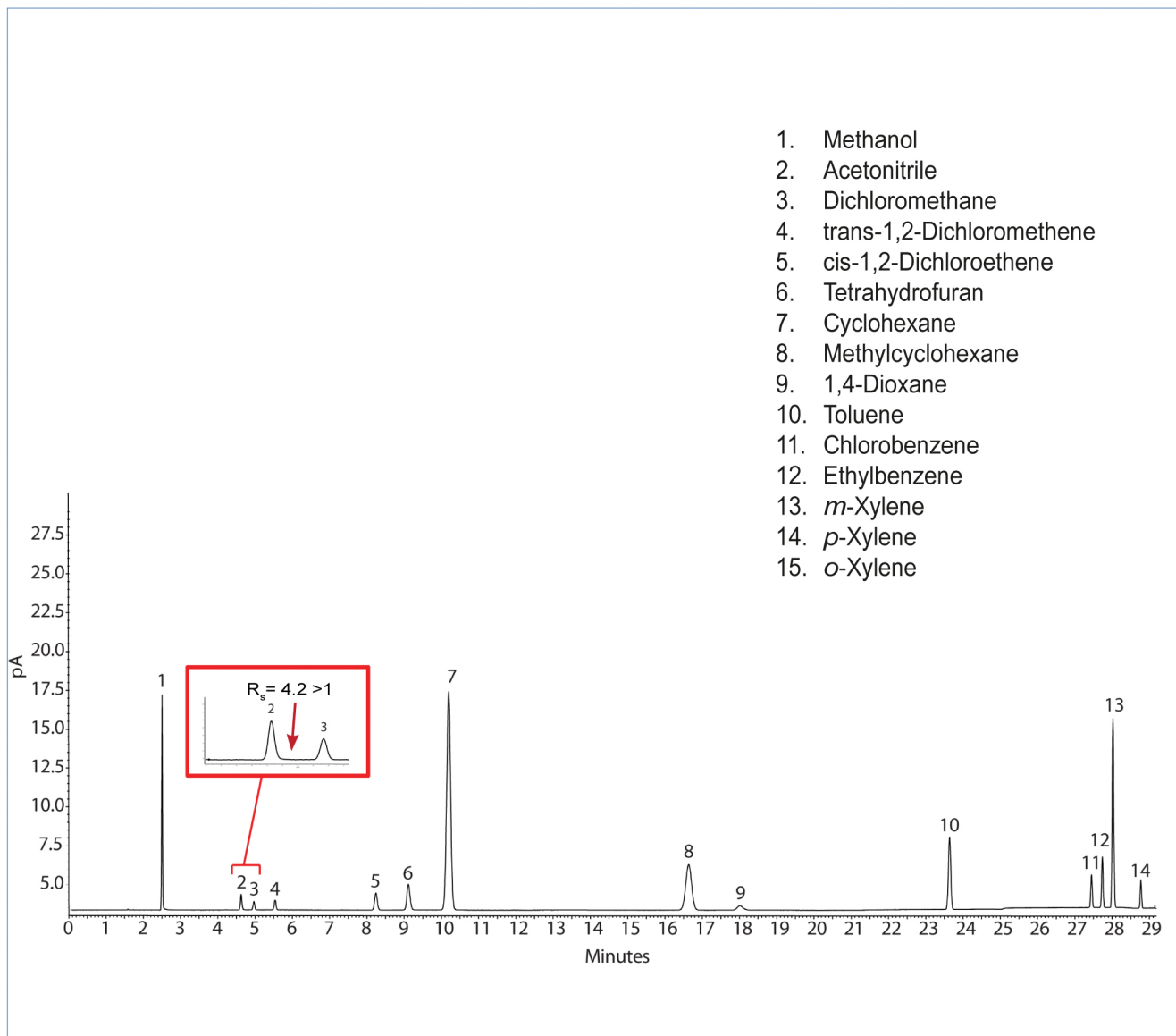

SAPIENS-624MS

www.teknokroma.es

Teknokroma Analítica S.A.- Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es