

SAPIENS-624MS: USP Glycerin monograph; Related compounds System Suitability Solution

Column: SAPIENS-624MS: P/N TR-593035


Dimensions: 30 m x 0,53 mm x 3 µm

Injection: 0.5 µl System suitability Solution (0.5mg/ml) 220°C, split: 1:10

Carrier gas: He, 38 cm/s

Oven temperature: 100°C to 220°C (4 min) @ 7.5°C/min

Detector: FID @ 250°C

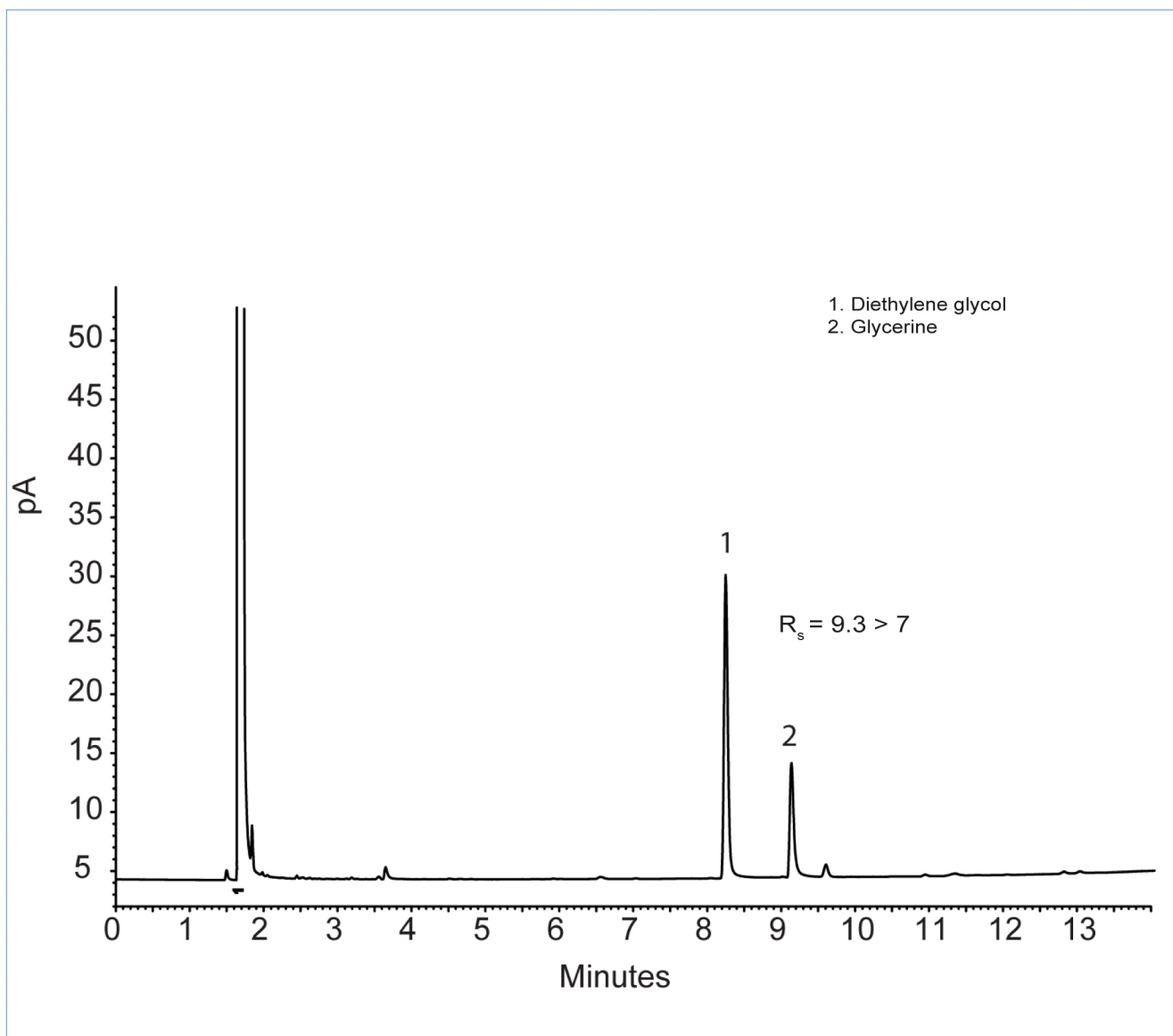

SAPIENS-624MS TKG 1306

www.teknokroma.esTeknokroma Analítica S.A.- Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es