


HYDROGENATED SOYA OIL

Column: **TR-CN100**, P/N TR-882162
Dimensions: 60m x 0.25mm x 0.20µm
Injection: 0.4µL Hydrogenated Soya Oil, split 1:50, 280°C
Carrier gas: H₂, 25 psi (172.2 KPa)
Oven program: 150°C(15min) @ 1.5°C/min to 220°C(5min)
Detector: FID, 280°C

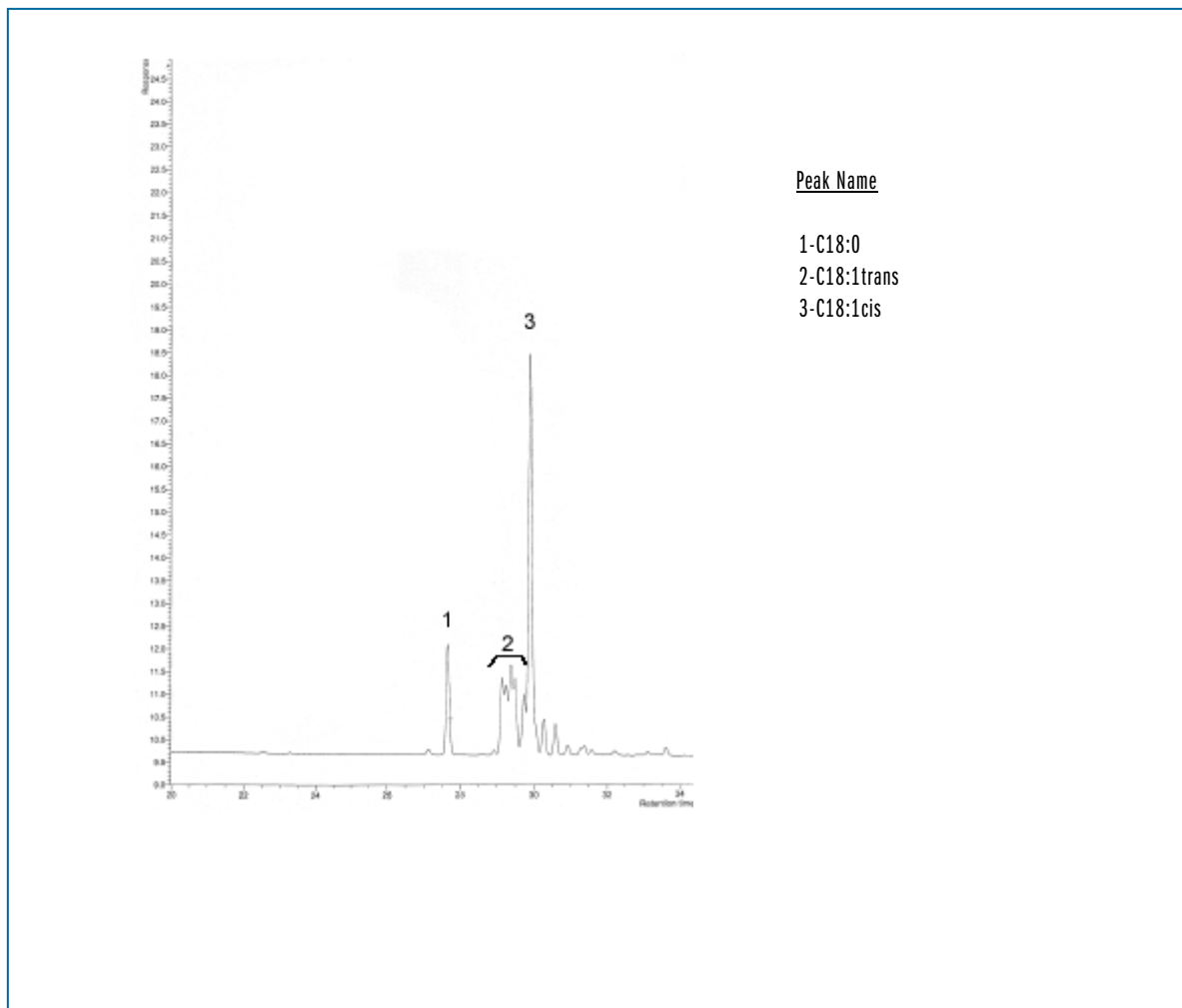

TKG 1024

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es