


LiteTouch® Fittings

Our LiteTouch Fittings are specifically designed to effectively grip tubing at two compression points, holding to high pressures with Fingertight convenience. It also prevents polymer tubing from twisting, a potential problem when using standard Fingertight fittings. LiteTouch Fittings are available for use with 1/32", 1/16", or 1/8" OD tubing sizes, and for 10-32 or 1/4-28 coned ports.

To avoid collapsing the ID of your tubing, the LiteTouch system can be used on hard tubing only, such as stainless steel and PEEK polymer tubing.

The LiteTouch Ferrule System is not recommended for repeated use in plastic ports.

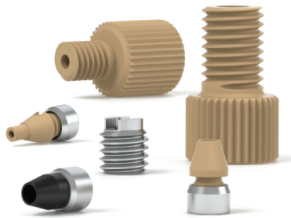

PEEK ZDV Unions

Cat.No.	Description
UP-F-354	LiteTouch® Stainless Steel FlushNut™, 10-32 Coned for 1/16" OD
UP-F-364	LiteTouch® Stainless Steel FlushNut™, 1/4-28 Coned, for 1/8" OD
UP-LT-100	LiteTouch® Ferrule w/Stainless Steel Lock Ring, 10-32 Coned, for 1/16" OD Natural - Single
UP-LT-100X	LiteTouch® Ferrule w/Stainless Steel Lock Ring, 10-32 Coned, for 1/16" OD Natural - 10 Pack
UP-LT-110	LiteTouch® PEEK Nut ,for 1/16" OD Single
UP-LT-110X	LiteTouch® PEEK Nut ,for 1/16" OD 10 Pack
UP-LT-132	LiteTouch® Micro Ferrule 10-32 Coned, for 1/32" OD
UP-LT-135	LiteTouch® SealTight™ Ferrule Black, 10-32 Coned for 1/16" OD Black - Single
UP-LT-200	LiteTouch® Ferrule 1/8" PEEK Natural
UP-LT-210	LiteTouch® Nut 1/8" PEEK Natural
UP-LT-215	LiteTouch® Short 1/8 Nut PEEK Natural Natural - Single
UP-LT-215B	LiteTouch® Short 1/8 Nut PEEK Natural Blue - Single
UP-LT-215BLK	LiteTouch® Short 1/8 Nut PEEK Natural Black - Single
UP-LT-215G	LiteTouch® Short 1/8 Nut PEEK Natural Green - Single
UP-LT-215GY	LiteTouch® Short 1/8 Nut PEEK Natural Gray - Single