

## TRB-1ms

Column: **TRB-1ms**, P/N TR-510262


Dimensions: 60m x 0.25mm x 0.25µm

Injection: 1µL Test MX5 (10 to 20 ng/comp. on column), split 1:100, 280°C

Carrier gas: H<sub>2</sub>, constant pressure 25 psi (172 KPa).

Oven temperature: 100°C @ 6°C/min to 325°C(5 min)

Detector: FID, 340°C

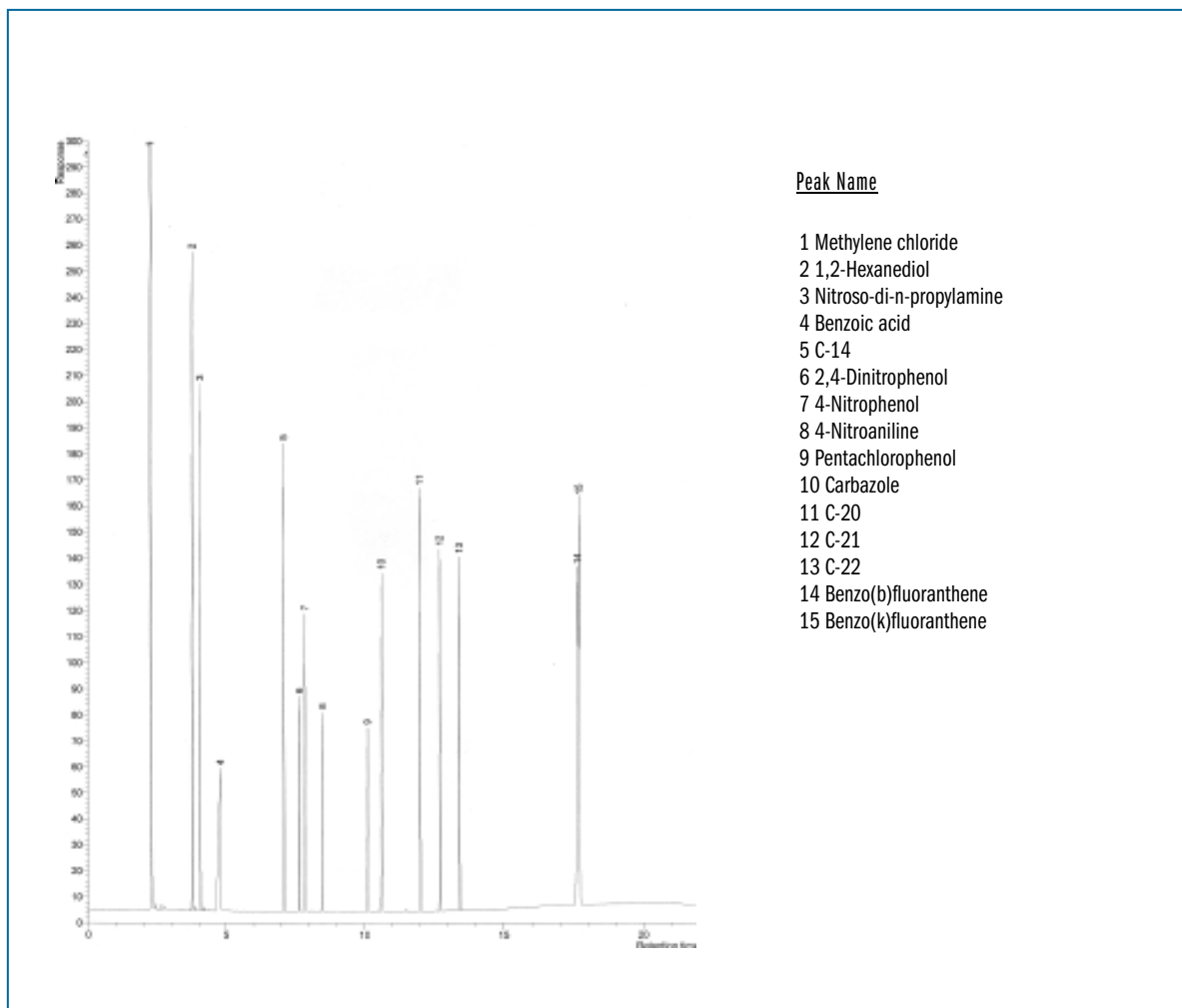

TKG 1104

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es