


## IRGANOX 1010

Column: **TBR-5ht**, 15m X 0,25 mm X 0,10 µm, P/N TR-620112  
Injection: 1µl (Irganox 1010, 12mg/ml chloroform), split (1:60), 370°C  
Carrier gas: H<sub>2</sub>, 6psi (41,3 kPa)  
Oven temp.: 150°C to 380°C (10 min.) @ 30°C/min.  
Detector: FID to 390°C

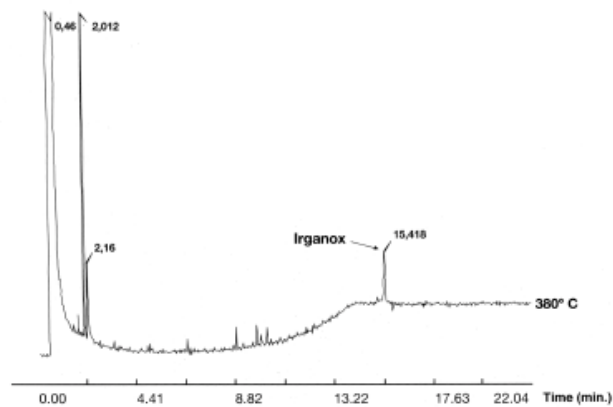

TKG 1110

[www.teknokroma.com](http://www.teknokroma.com)