

Sterols

Column: **TBR-Sterol**, 30m X 0,22 mm X 0,22 μ m, P/N TR-182238


Oven Temp.: 265°C

Injector: 280°C

Carrier gas: H₂, 18 psi (124 kPa)

Injection: 0,5 μ l sterols standard, (25 mg/ml.) split(1:100)

Detector: FID 300°C

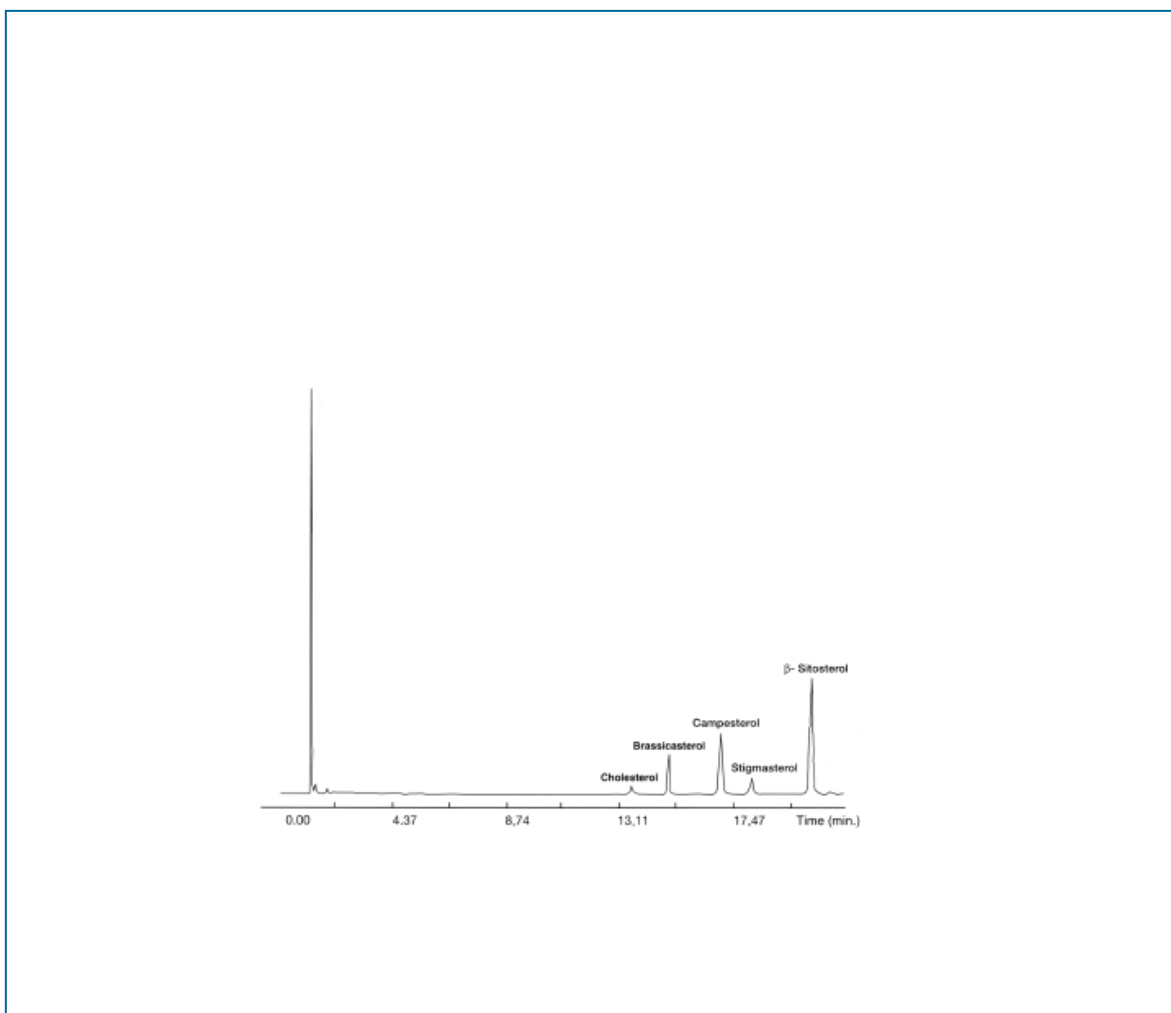

TKG 1111

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es