


DRUGS IN URINE

Column: **TRB-5ms**, P/N TR-520129
Dimensions: 25m x 0.20mm x 0.11µm
Injection: 250°C, 1 µl splitless (BSTFA Derivatives in ACN)
Carrier gas: He, 15 psi (103.3 kPa)
Oven temperature: 60°C (1') to 180°C (1') @ 10°C/min. to 220°C @ 10°C/min.
Detector: FID, 280°C

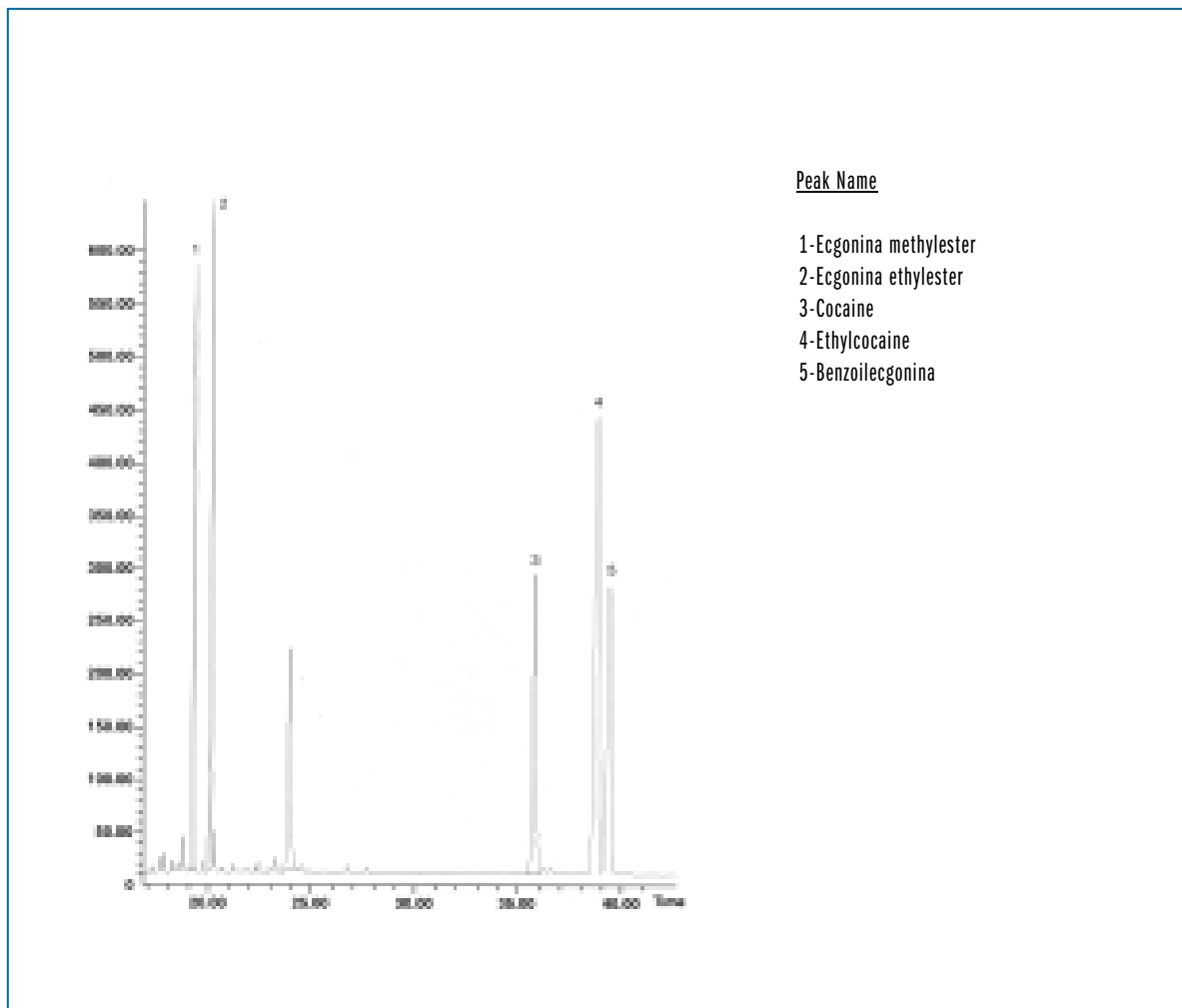

TKG 1067

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es