

**PUFA III**

Column: **SupraWax-280**, 15m x 0.10mm x 0.10µm (P/N: TR-830111)


Injection: 280°C, split 200:1, precission liner

Detector: FID, 280°C

Carrier Gas: H<sub>2</sub>, 45 psi

Oven: 100°C (0.5min) to 280°C @ 50°C/min (2min)

Sample: 0.1µL PUFA III – Partially Hydrogenate Menhaden Oil diluted to 100mg/ml in hexane

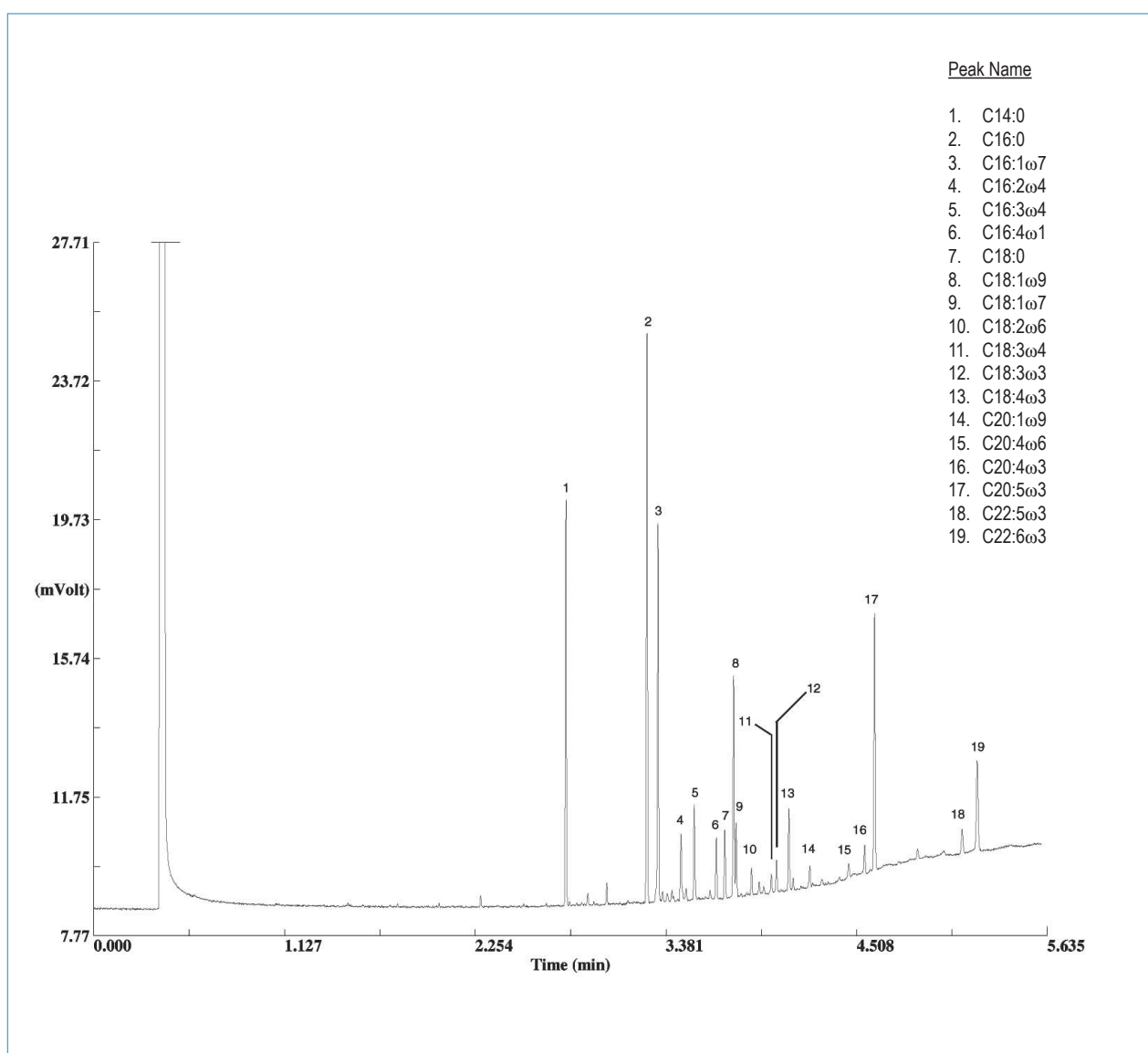

**TKG 1253**

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es