


SEPARATION OF FAMES

Column: **SupraWax-280**, P/N TR-230232
 Dimensions: 30m x 0.25mm x 0.25µm
 Injection: 1µL SP-47885-U (10 mg/mL), split 1:100, 260°C
 Carrier gas: He, 30 psi (206.7 KPa)
 Oven program: 50°C(2min) @ 4°C/min to 220°C(15min)
 Detector: FID, 260°C

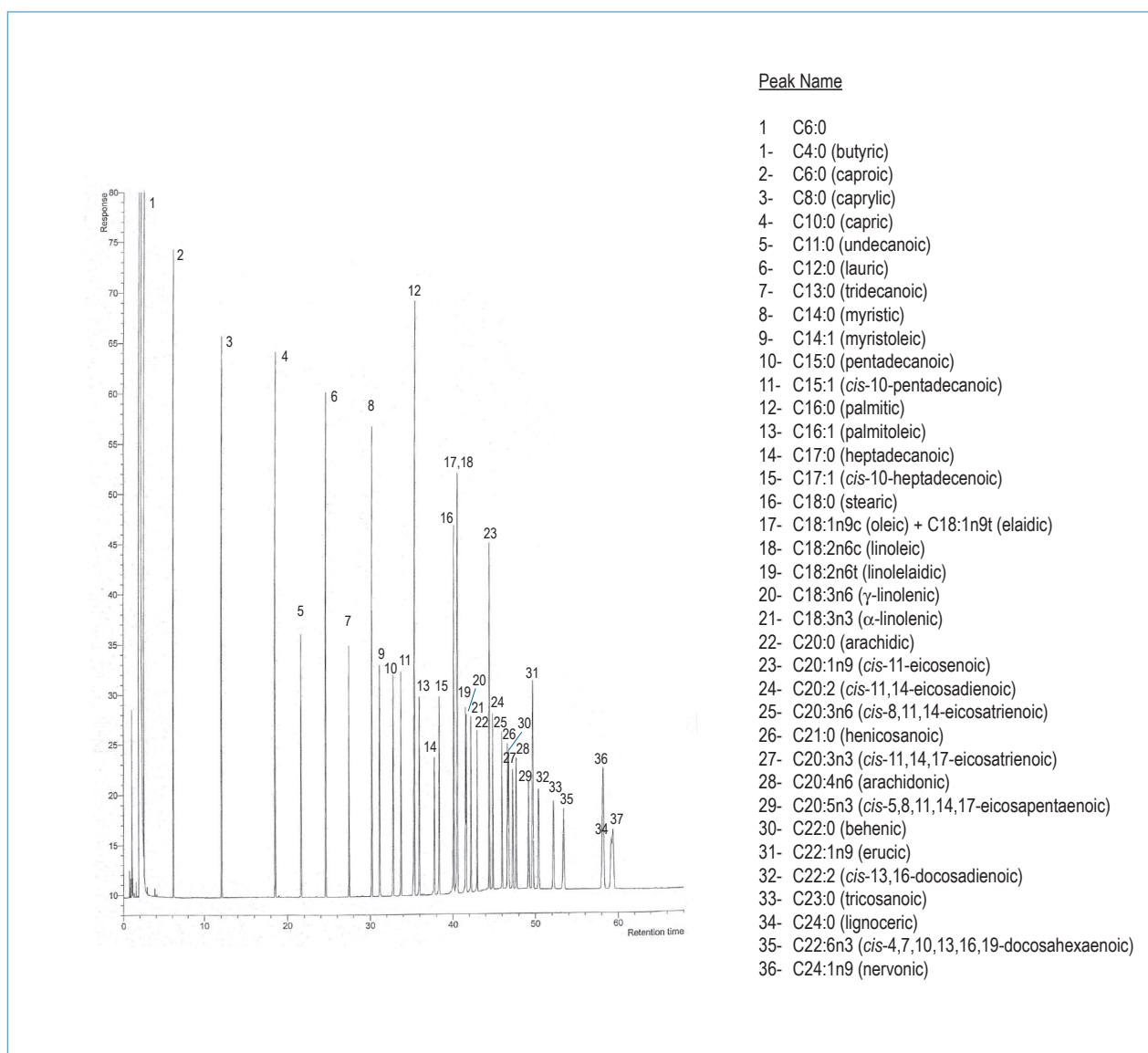

TKG 1237

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es