

GLYCOLS IN WINE

Column: **SupraWAX-280**, P/N TR-831035

Size: 30m x 0.53mm x 1.0 μ m


Injection: 1 μ L split 1:50, 220°C

Carrier gas: He, 5.4 psi

Program temperature: 100°C (3min) @ 4.5°C/min to 220°C (10min)

Detector: FID, 260°C

Chromatogram provided by R. Franquet and Joan Garcia from INCAVI (Vilafranca del Penedés, Barcelona)

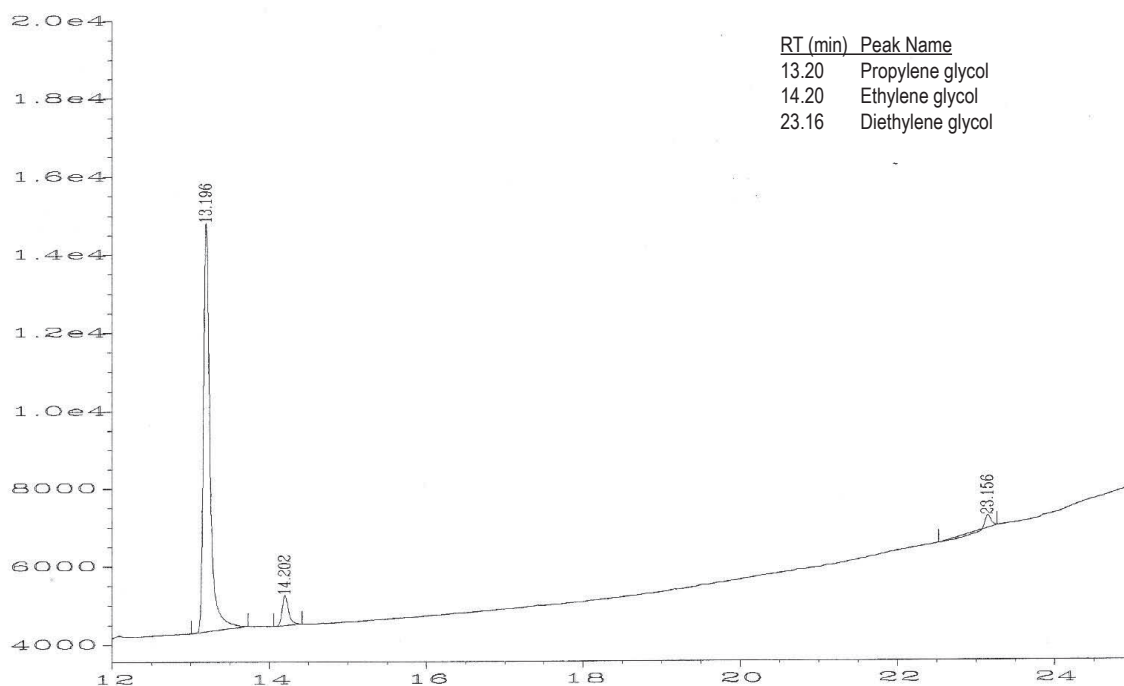

TKG 1241

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es