


## SEPARATION OF FAMES

Column: **TR-CN100**, P/N TR-882162  
 Dimensions: 60m x 0.25mm x 0.20µm  
 Injection: 1µL SP-47885 (10 mg/mL), split 1:100, 260°C  
 Carrier gas: He, 20.15 cm/s at 185°C  
 Oven program: 90°C(7min) @ 5°C/min to 240°C(15min)  
 Detector: FID, 260°C

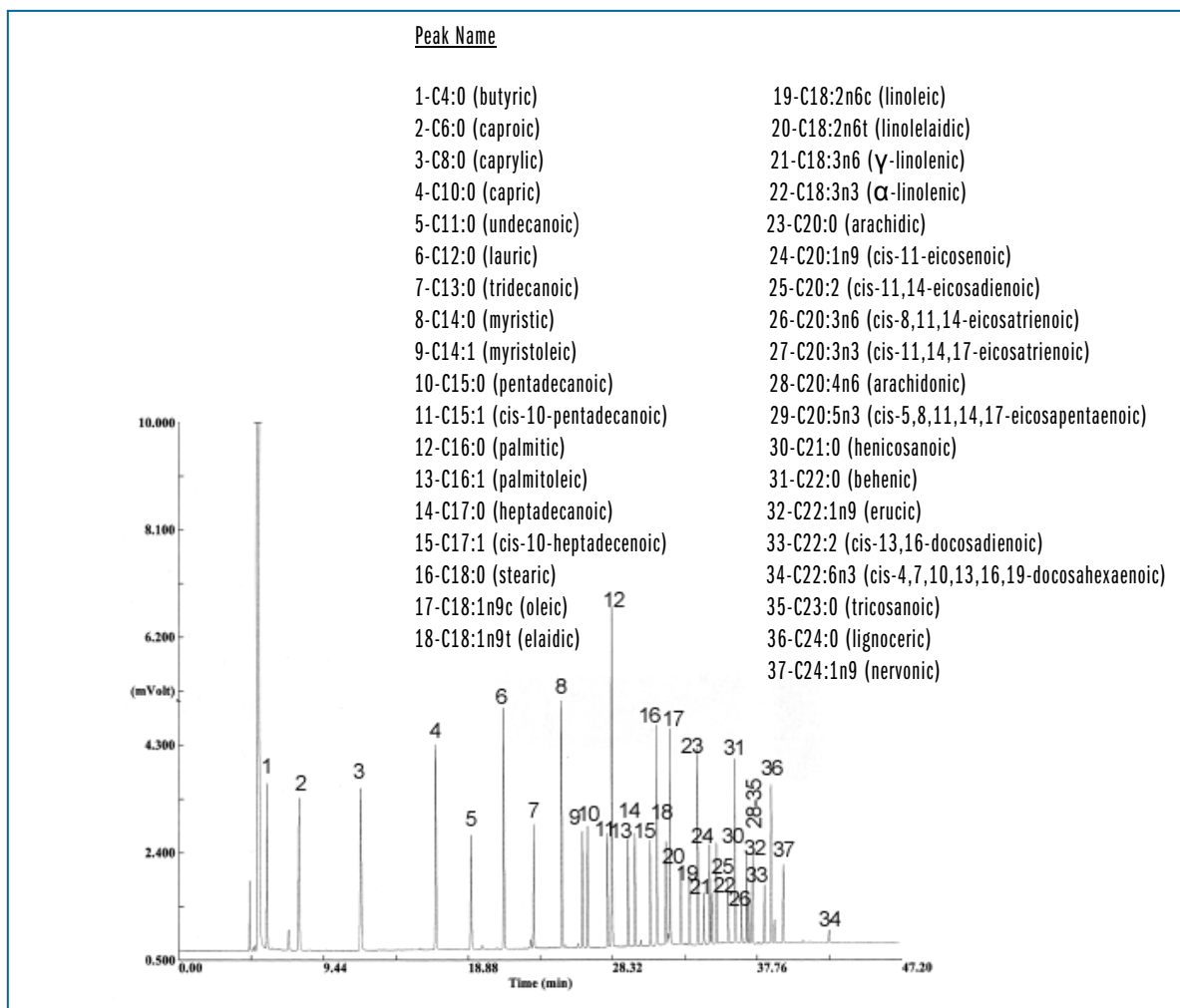

TKG 1025

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es