


ANALYSIS OF CIS-TRANS ISOMERS

Column: **TR-CN100**, P/N TR-882133
Dimensions: 30m x 0.32mm x 0.20mm
Injection: 1 µL isomers standard, split
Carrier gas: H₂, 4.5 psi (31 KPa)
Oven temperature: 140°C @ 4°C/min to 190°C
Detector: FID, 250°C

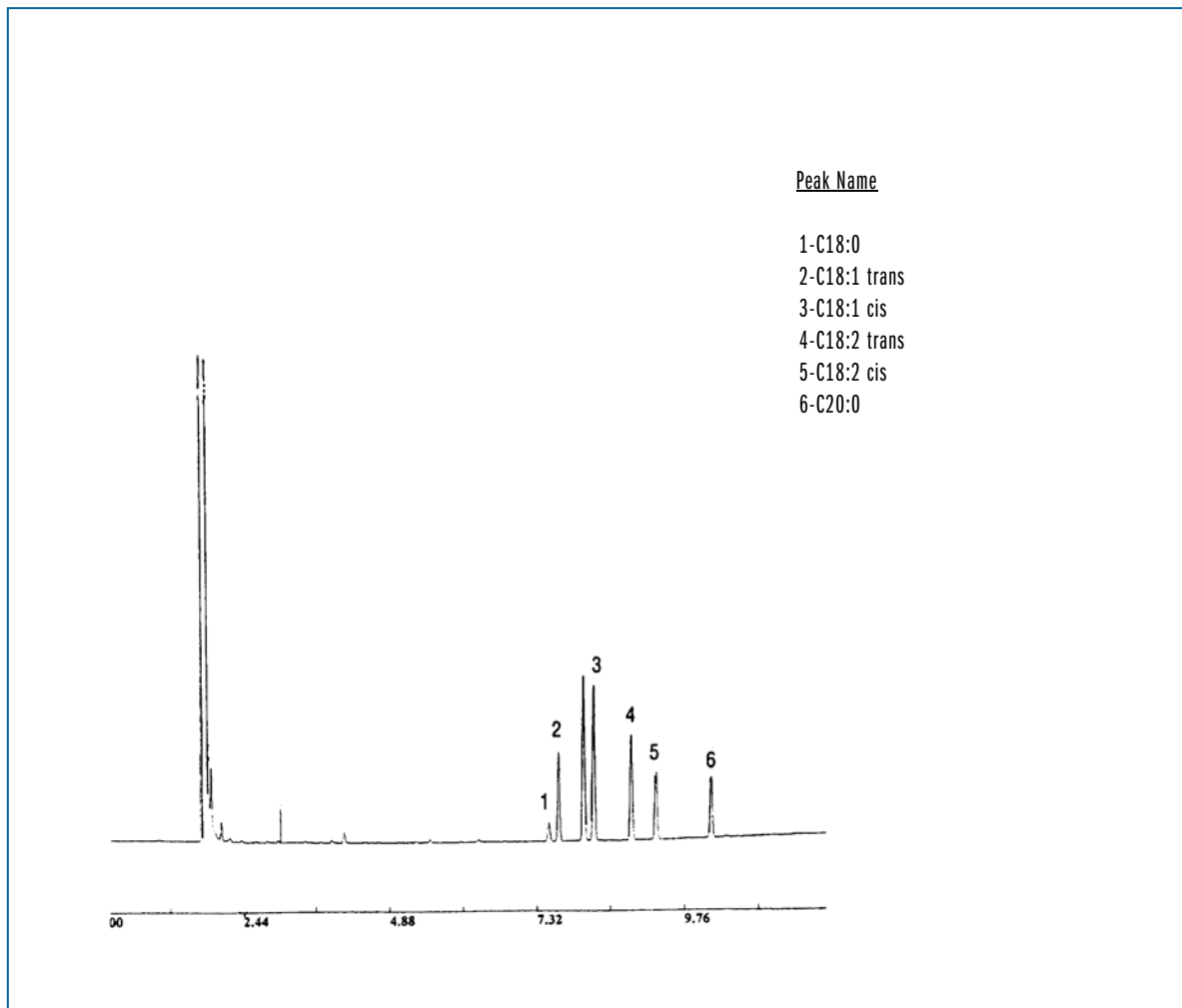

TKG 1185

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es