


## ANALYSIS OF METHYL ESTERS

Column: TR-CN100, P/N TR-882113  
 Dimensions: 15m x 0.32mm x 0.20µm  
 Injection: 1 µL FAMES standard, split  
 Carrier gas: H<sub>2</sub>, 2.4 psi (16.6 KPa)  
 Oven temperature: 140°C @ 3°C/min to 185°C  
 Detector: FID, 250°C

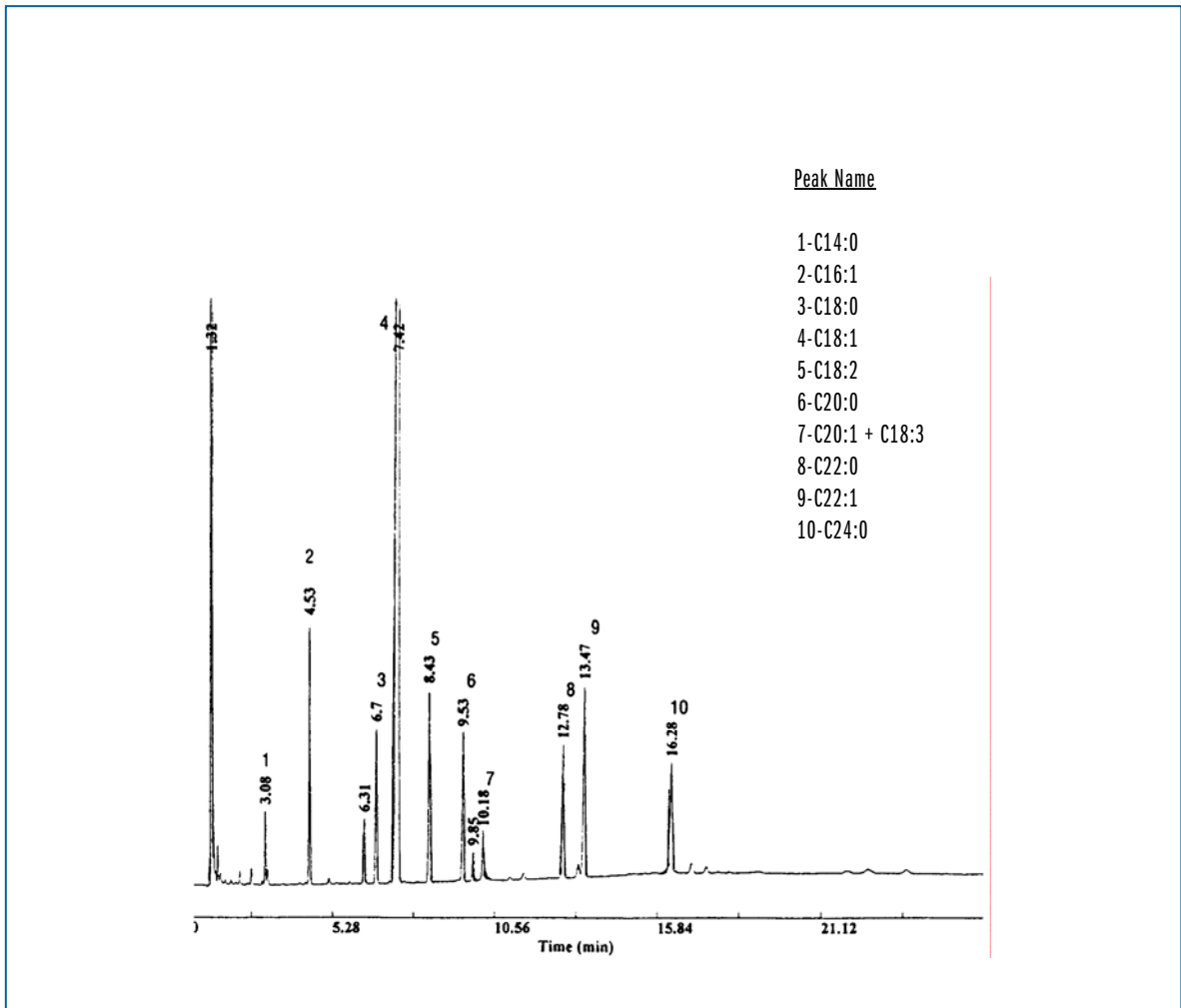

TKG 1187

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es