

MAXIMUM SEPARATION OF CIS-TRANS FAMES

Column: **TR-CN100**, P/N TR-882192


Dimensions: 100m x 0.25mm x 0.20µm

Injection: 0.5µL cis/trans FAMES standard (10mg/mL), split 1:80, 260°C

Carrier gas: H₂, constant flow 0.8mL/min

Oven Temperature: 175°C

Detector: FID, 260°C

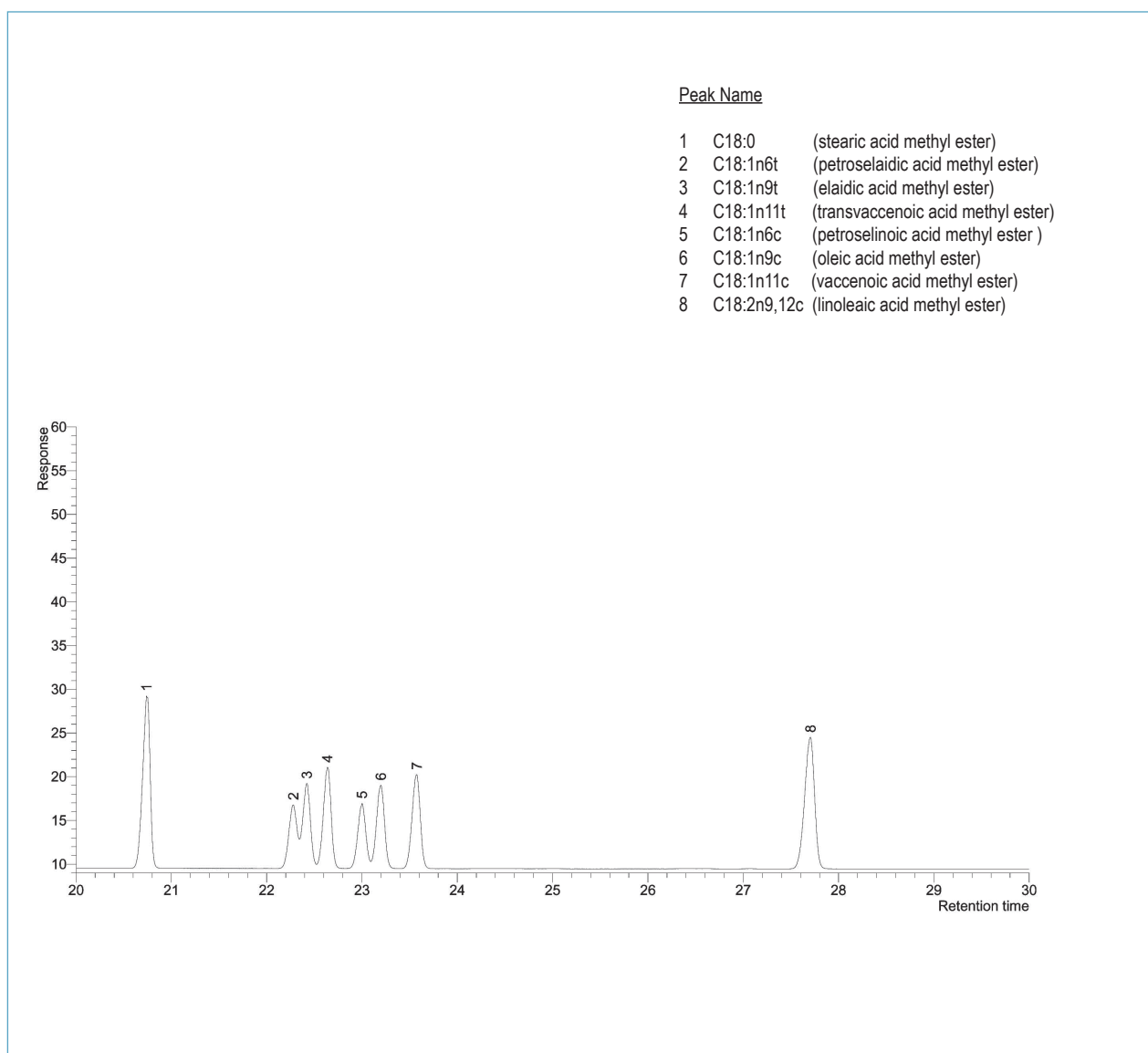

TKG 1255

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es