

## IMPURITIES IN CYCLOHEXANE

Column: **TRB-1**, P/N TR-111226


Dimensions: 25m x 0.15mm x 1.2  $\mu$ m

Injection: 1 $\mu$ L Cyclohexane, split 1:100, 260°C

Carrier gas: H<sub>2</sub>, constant pressure 22 psi (0.8 mL/min)

Oven temperature: 60°C(6 min) @ 20°C/min to 160°C(2 min)

Detector: FID, 300°C

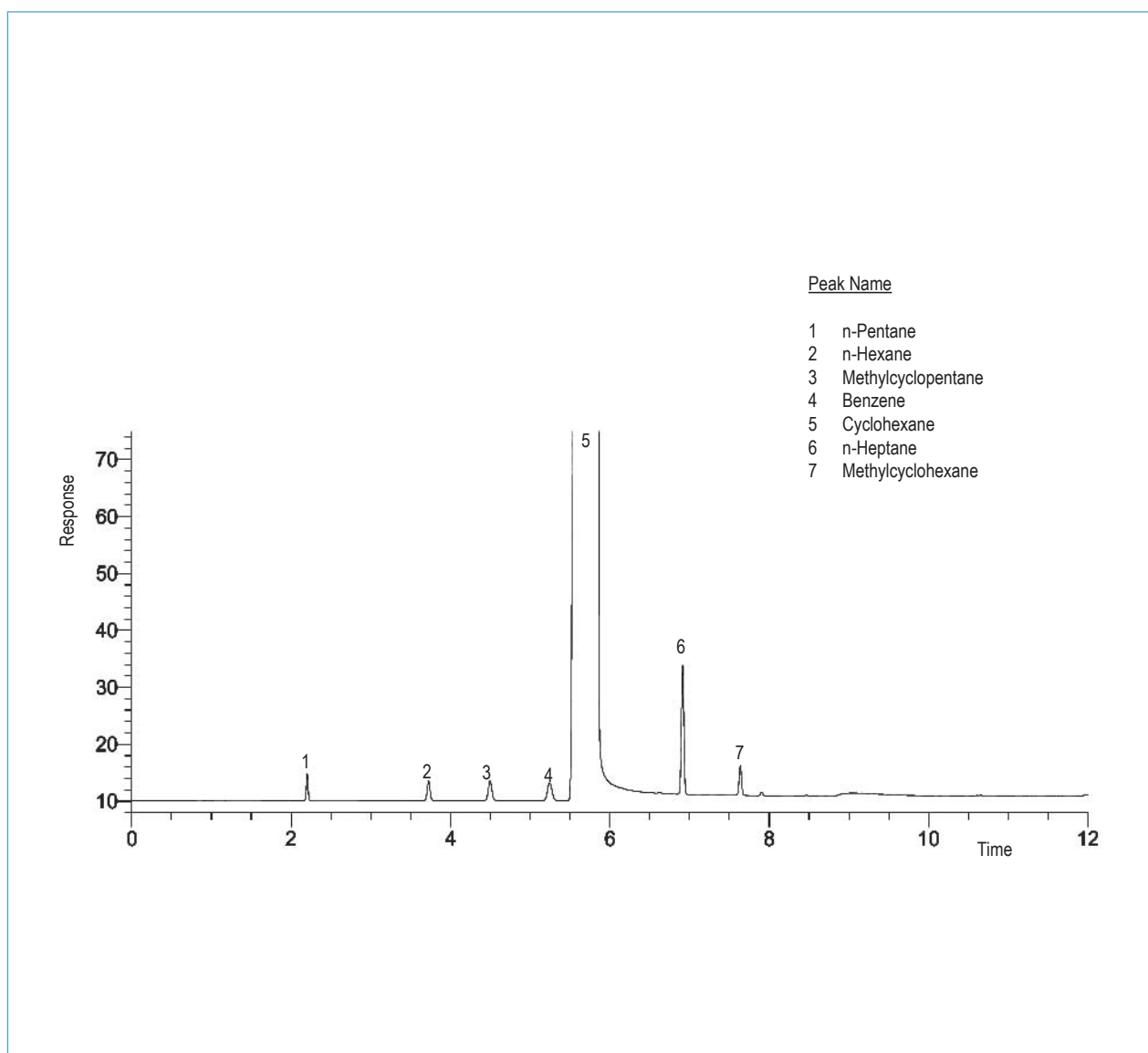

TKG 1195

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: [comercial@teknokroma.es](mailto:comercial@teknokroma.es)