

VOLATIL SOLVENTS NO HALOGENATED

Column: **TRB-1**, P/N TR-111062


Size: 60m x 0.25mm x 1.0 μ m

Injection: 1 mL Headspace (70°C, 20min), split 1:5, 250°C

Carrier Gas: He, constant flow 1mL/min

Program Temperature: 40°C(2 min) @ 8°C/min to 240°C(10 min)

Detector: MS, ion source 200°C, Interfase 250°C, scan 20-400

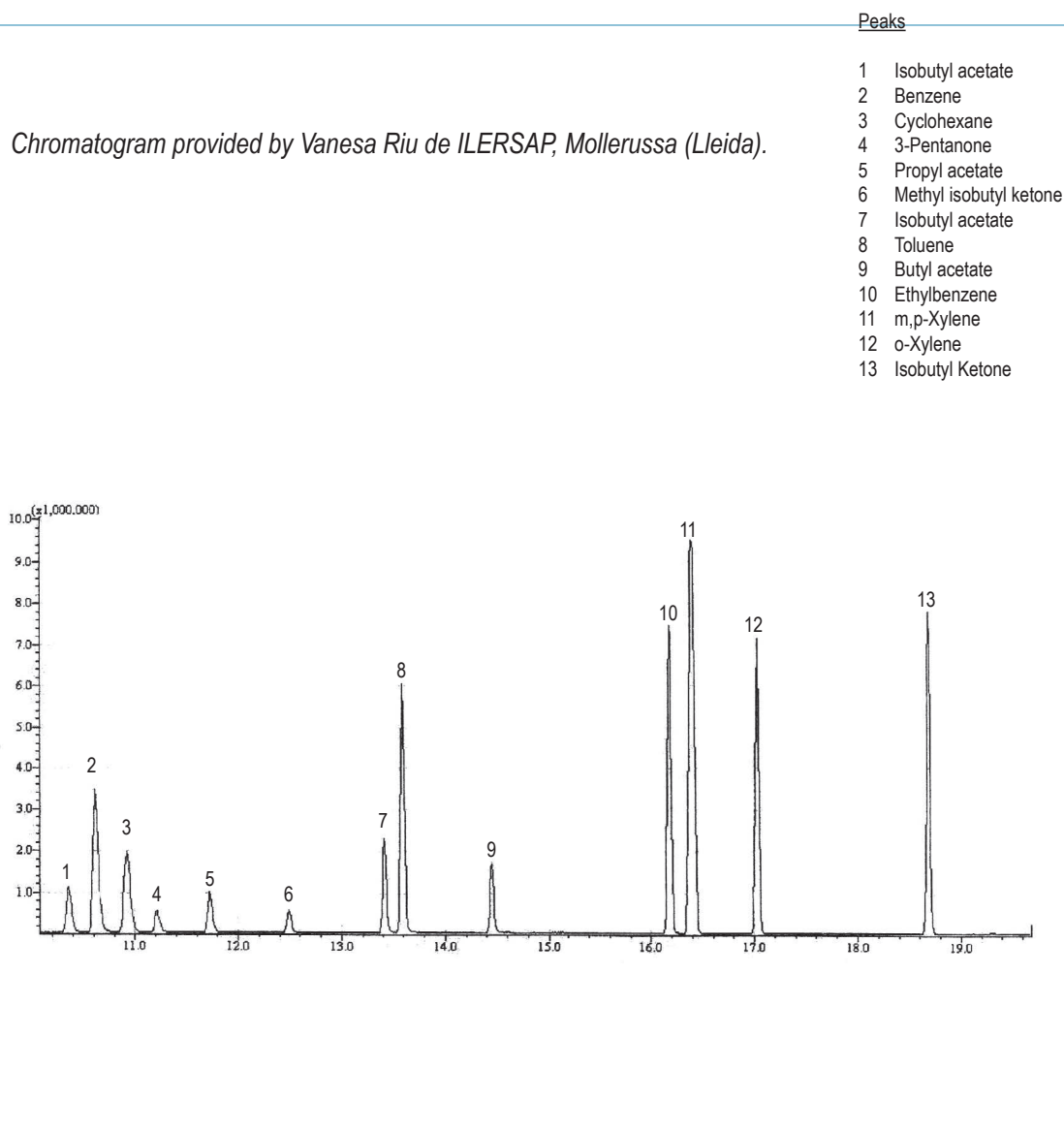

TKG 1220

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es