

COMMON INDUSTRIAL SOLVENTS

Column: **TRB-1**, P/N TR-111033


Size: 30m x 0.32mm x 1.0 μ m

Injection: 0.01 μ L Neat solvents, split 1:300, 200°C

Carrier Gas: H₂, 7 psi

Program temperature: 30°C @ 8°C/min to 140°C(2min)

Detector: FID, 200°C

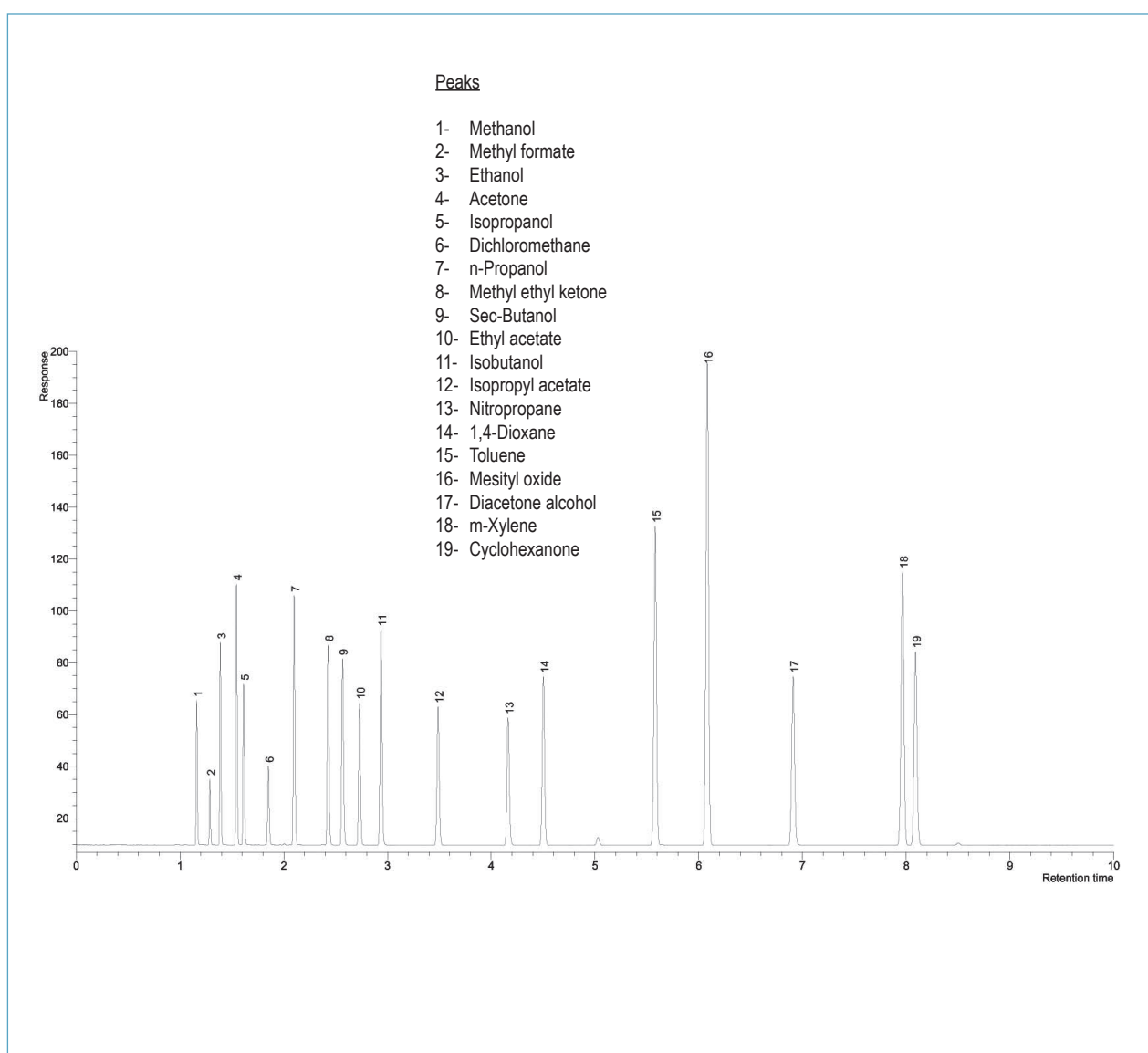

TKG 1234

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es