

## PESTICIDES

Column: **TRB-5**, P/N TR-120469


Dimensions: 60m x 0.20mm x 0.40µm

Injection: 1µL standard Organochloride Pesticide Mix AB#2, splitless(1min), 270°C

Carrier gas: He, constant pressure 39 psi (268.7 Kpa)

Oven temperature: 90°C(3min) @ 30°C/min to 215°C(40min) @ 5°C/min to 275°C(30min)

Detector: ECD, 300°C

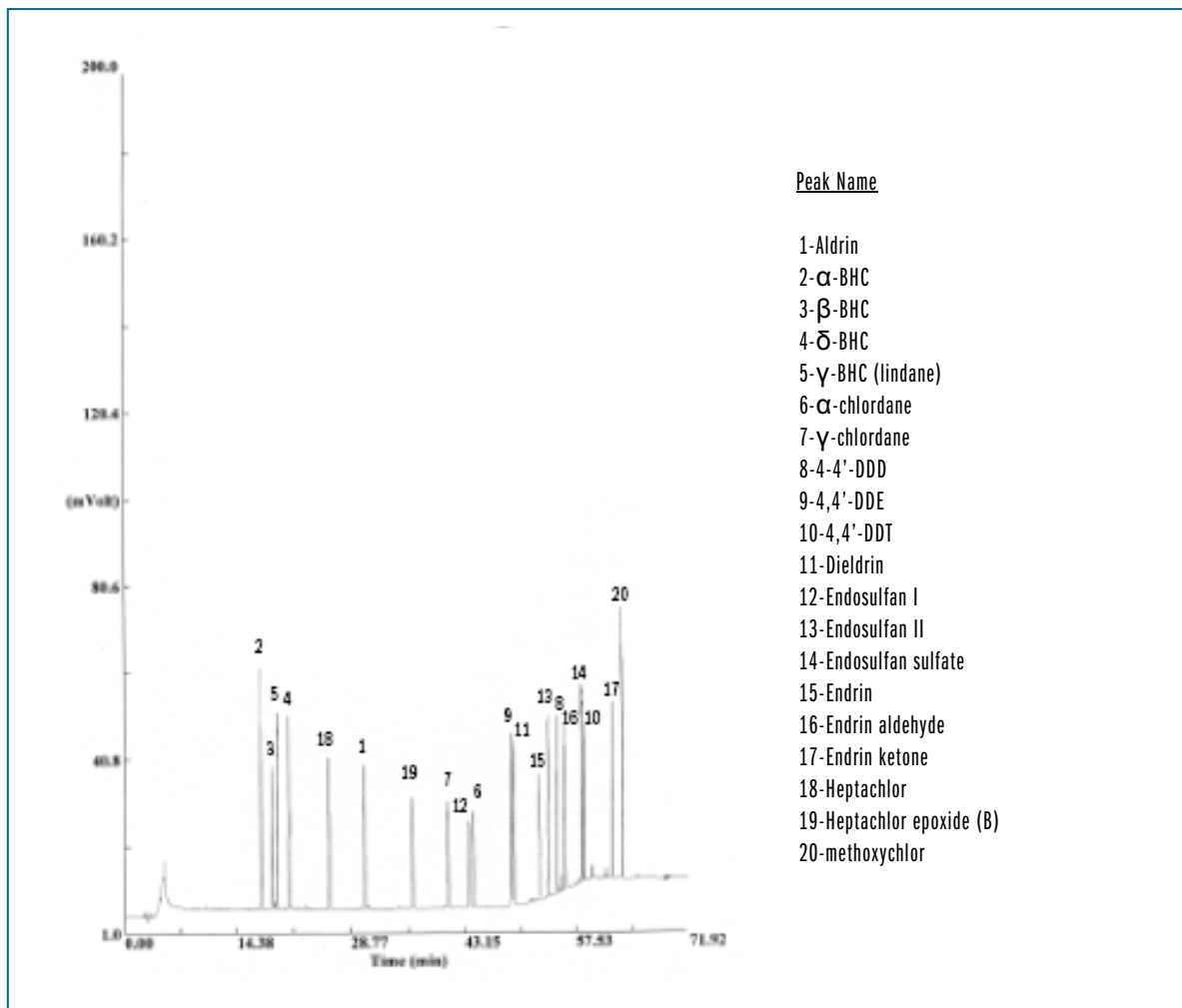

TKG 1055

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es