


## ANALYSIS OF NONYLPHENOLS

Column: **TRB-5**, P/N TR-120262  
Dimensions: 30m x 0.25mm x 0.25µm  
Injection: 1µL, split  
Carrier gas: H<sub>2</sub>, 50 cm/s (110°C)  
Oven temperature: 110°C @ 20°C/min to 220°C(1min) @ 4°C/min to 300°C  
Detector: FID, 310°C

*Chromatogram provided by Dr. Caixach from Laboratori Espectrometria Masses, CSIC, Barcelona.*

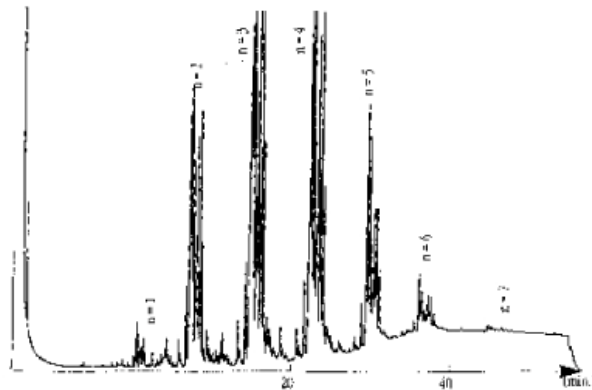

TKG 1151

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es