


## ANALYSIS OF AROCLORS

Column: **TRB-5**, P/N TR-122168  
Dimensions: 60m x 0.22mm x 0.20µm  
Injection: splitless  
Carrier gas: H<sub>2</sub>, 150 Kpa  
Oven temperature: 80°C(3,1min) @ 50°C/min to 190°C(5min)  
@ 1°C/min to 230°C(4min) @ 4°C/min to 260°C  
Detector: ECD, 350°C

*Chromatogram provided by C. Rodriguez and L. Comella from Institut Químic de Sarrià, Barcelona.*

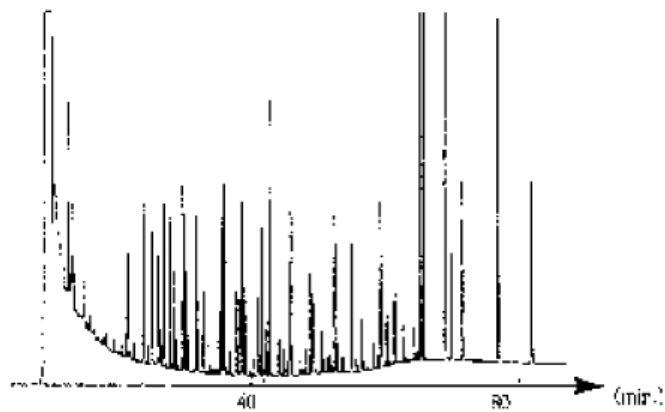

TKG 1152

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es