


PHENOLS EPA 604

Column: **TRB-5**, P/N TR-120232
Dimensions: 30m x 0.25mm x 0.25µm
Injection: 1 µL, split, 2 to 6 ng/comp, 250°C
Carrier gas: H₂, 12 psi (82.68 KPa)
Oven temperature: 80°C(4min) @ 8°C/min to 250°C
Detector: FID, 280°C

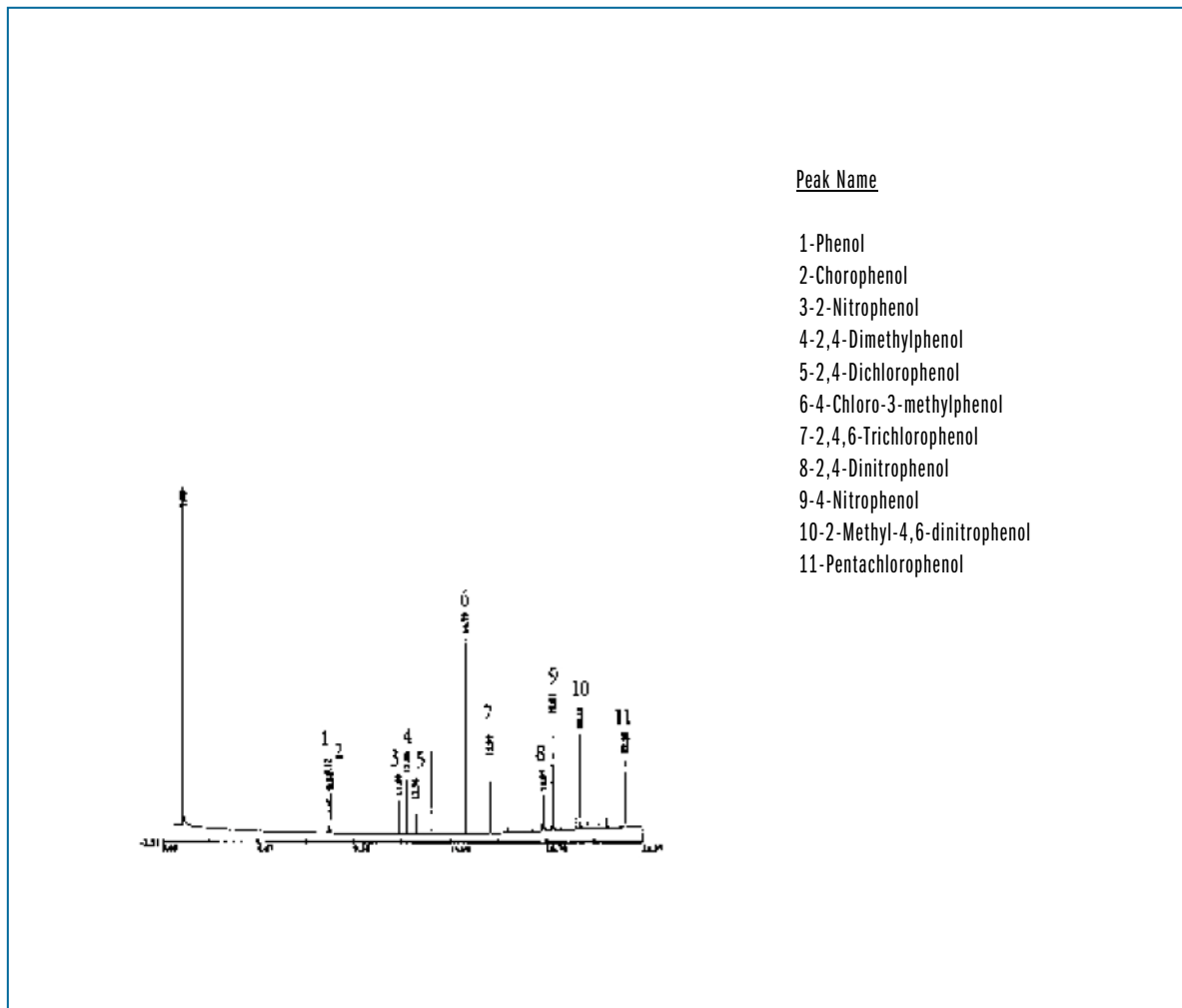

TKG 1153

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es