


ANALYSIS OF PESTICIDES

Column: **TRB-5**, P/N TR-120232
Dimensions: 30m x 0.25mm x 0.25 μ m
Injection: 1 μ L, split
Carrier gas: He, 14 psi (96.46 KPa)
Oven temperature: 150°C @ 5°C/min to 265°C
Detector: FID, 325°C

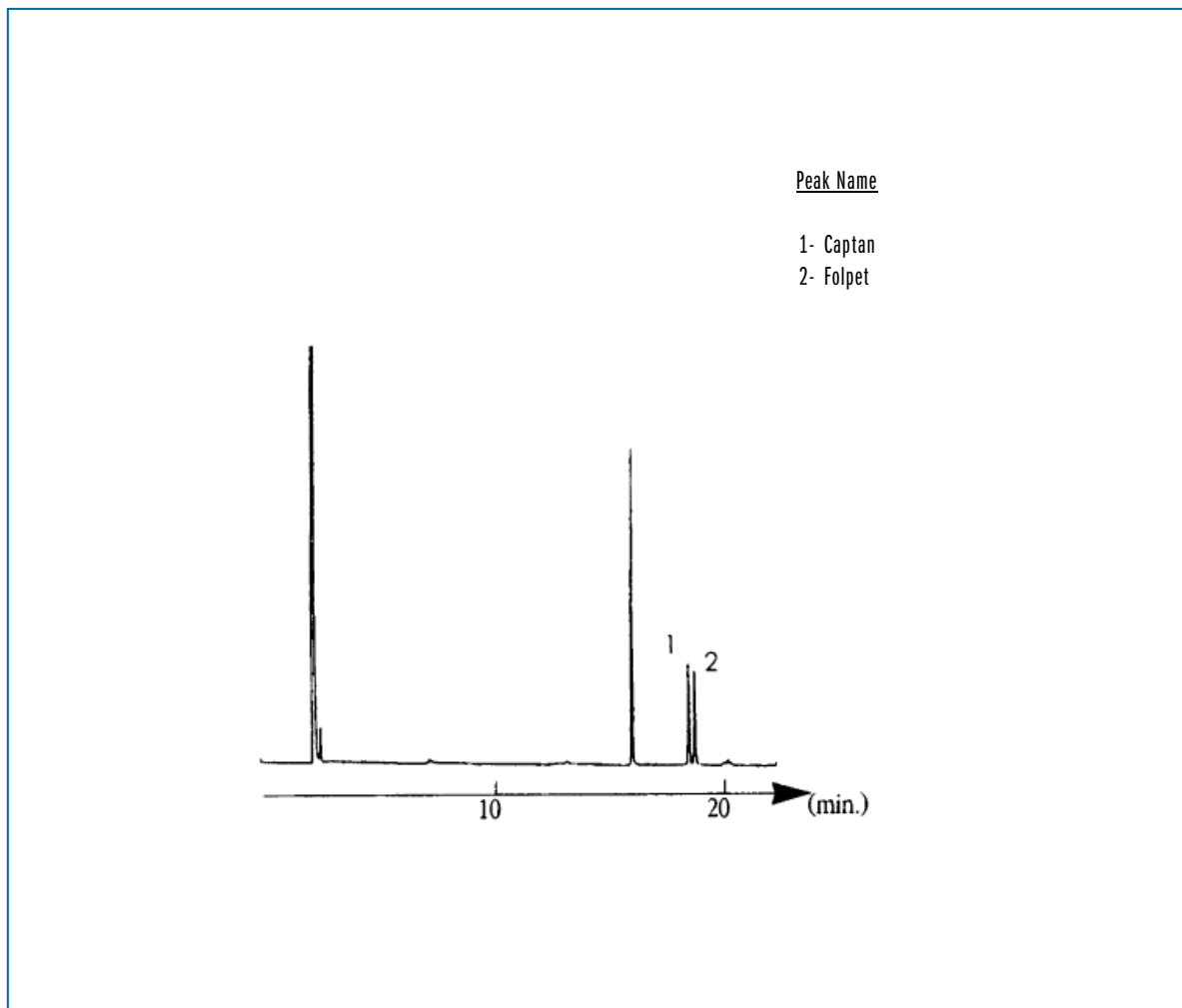

TKG 1155

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es