


## ANALYSIS OF HYDROCARBONS (GASOIL)

Column: **TRB-5**, P/N TR-120232  
Dimensions: 30m x 0.25mm x 0.25µm  
Injection: 1 µL, splitless 0.7min  
Carrier gas: H<sub>2</sub>, 50 cm/s (110°C)  
Oven temperature: 60°C(3 min) @ 4°C/min to 300°C  
Detector: FID, 305°C

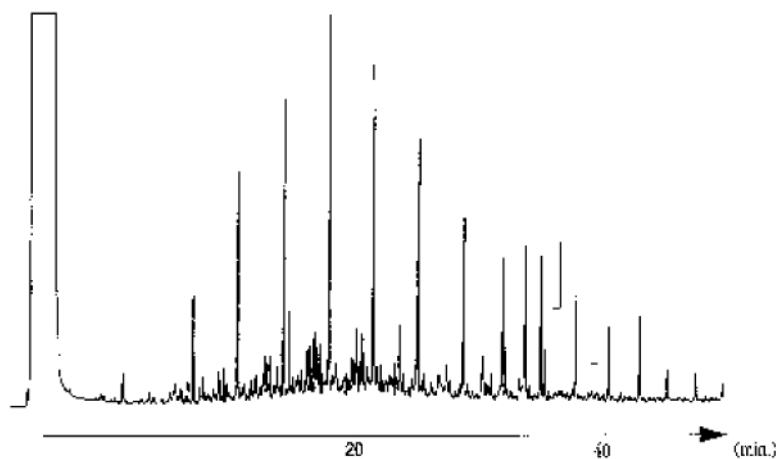

*Chromatogram provided by Dr. Caixach from Laboratori Espectrometria de Masses, CSIC, Barcelona*

TKG 1173

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es