

AMANTADINE HYDROCHLORIDE IMPURITIES

Column: **TRB-5 AMINE**, P/N TR-211035

Dimensions: 30m x 0.53mm x 1.0m


Injection: 2 µl (split 1:50), 220°C

Carrier gas: He, 4.2 psi (28.9 kPa)

Oven temperature: 70°C (5') to 250°C (20min) @ 10°C/min.

Detector: FID, 300°C

Sample: Test solution according to USP 25

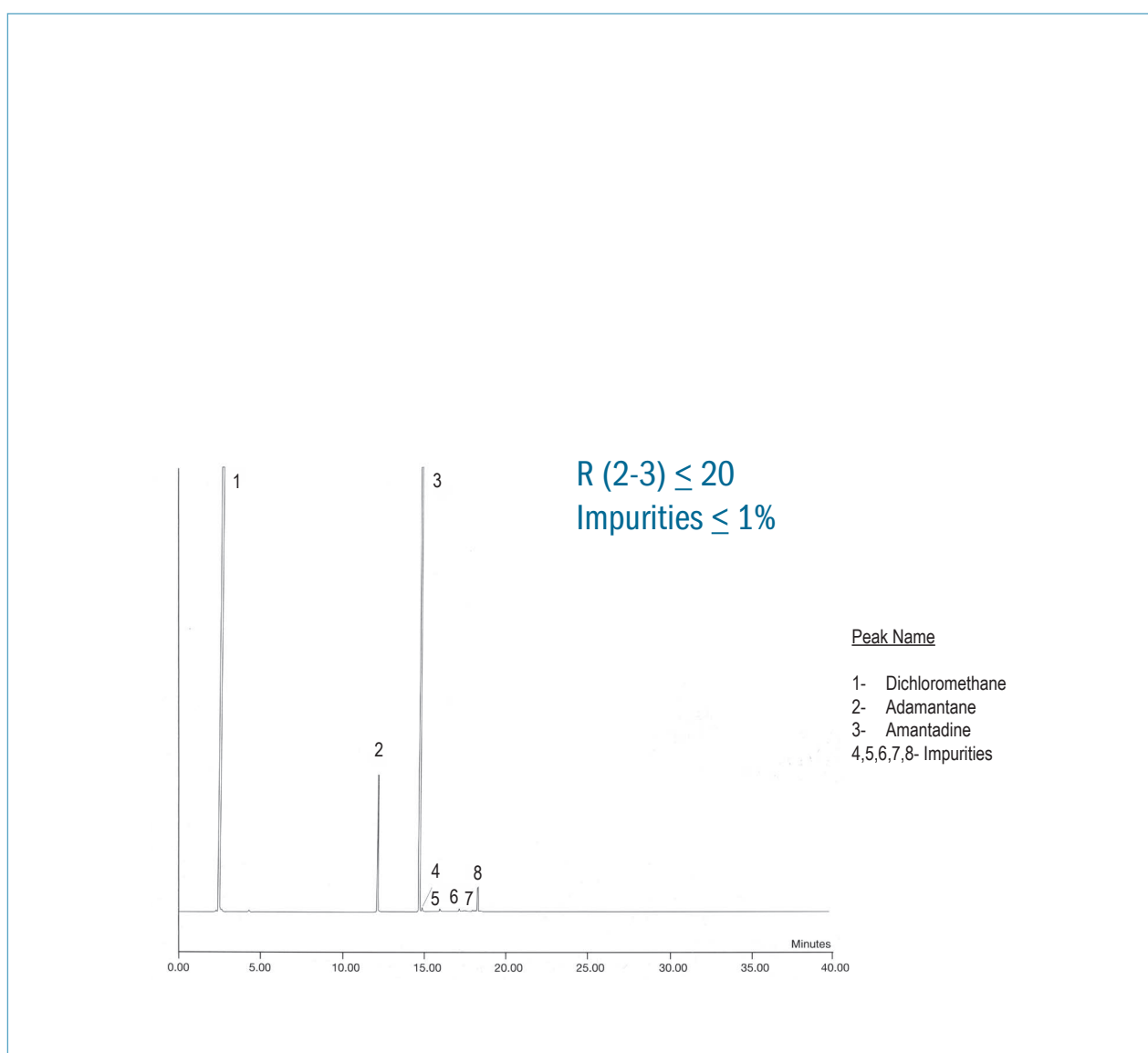

TKG 1074

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es