

AMINES

Column: TRB-5A, P/N TR-213035

Size: 30m x 0.53mm x 3.0 μ m


Injection: 1 μ L standard (wet needle), split 1:100 (50 ng/comp), 280 $^{\circ}$ C

Carrier Gas: H₂, constant pressure, 3 psi

Program temperature: 35 $^{\circ}$ C(3min) @ 10 $^{\circ}$ C/min to 225 $^{\circ}$ C(2min)

Detector: FID, 300 $^{\circ}$ C

RT (min)	Peak Name
1.58	Ethyl amine
2.0	Isopropylamine
4.25	Isobutylamine
5.14	n-Butylamine
14.63	n-Octylamine
16.49	n-Nonylamine
18.22	n-Decylamine

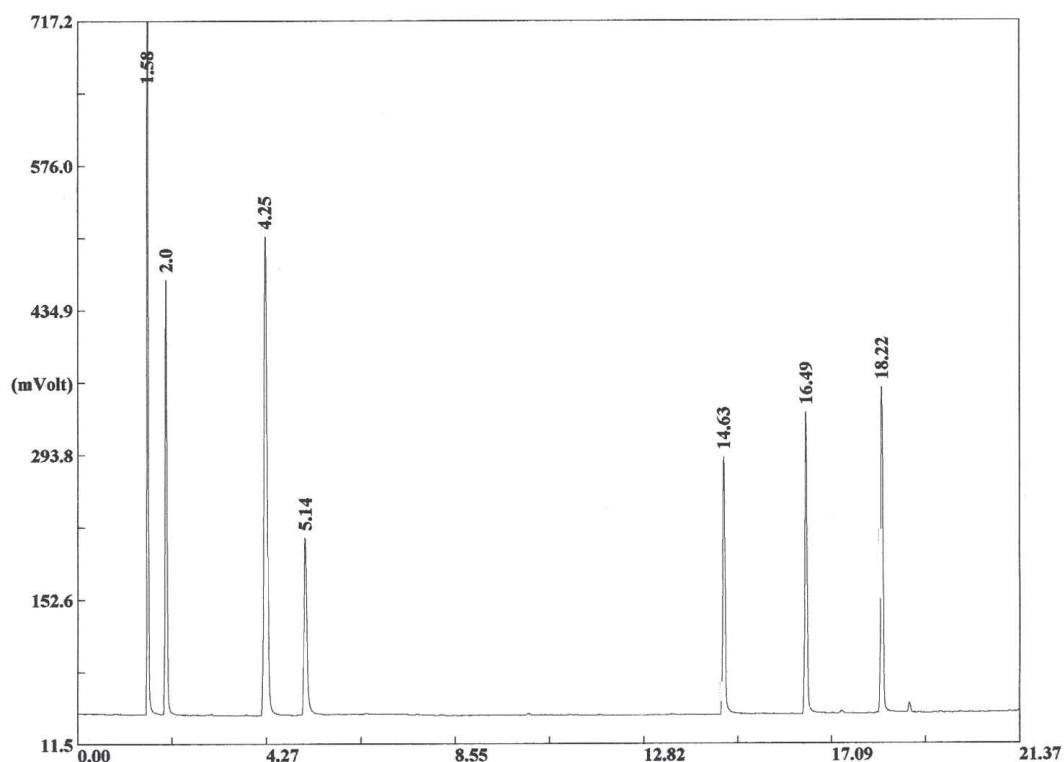

TKG 1212

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es