


Aromatic Amines (Azo dye Compounds)

Column: TRB-5A, P/N: TR-210532
 Dimensions: 30m x 0.25mm x 0.50um
 Injection: 260°C, split 1:15
 Carrier Gas: He, constant flow @ 1ml/min
 Oven temperature: 100°C (2min) to 310°C (2min) @ 15°C/min
 Detector: MS (SIM)
 Interface temperature: 280°C
 Ionization mode: EI Sample: 1ul Mix Aromatic Amines + Aniline + 1,4-Phenylenediamine (10ppm each compound)

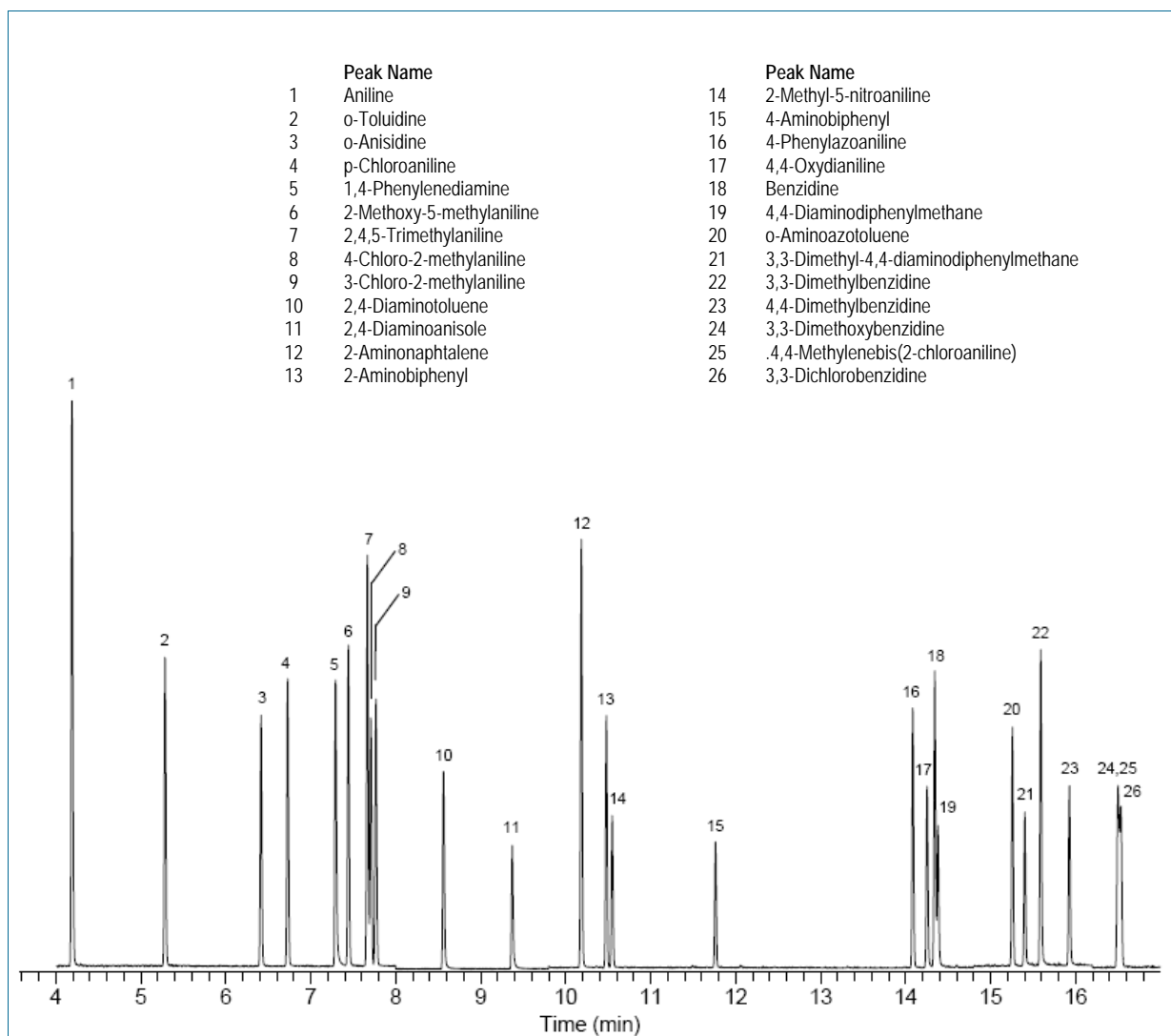

TKG 1279