

## Aromatic Amines (Azo dye Compounds)

Column: **TRB-5A**, P/N: TR-210232

Dimensions: 30m x 0.25mm x 0.25µm

Injection: 260°C, split 1:15

Carrier Gas: He, constant flow @ 1ml/min


Oven temperature: 55°C to 150°C (2min) @ 30°C/min to 280°C (2min) @ 8°C/min

Detector: MS (SIM)

Interface temperature: 280°C

Ionization mode: EI

Sample: 1µl Mix Aromatic Amines + Aniline + 1,4-Phenylenediamine (10ppm each compound)

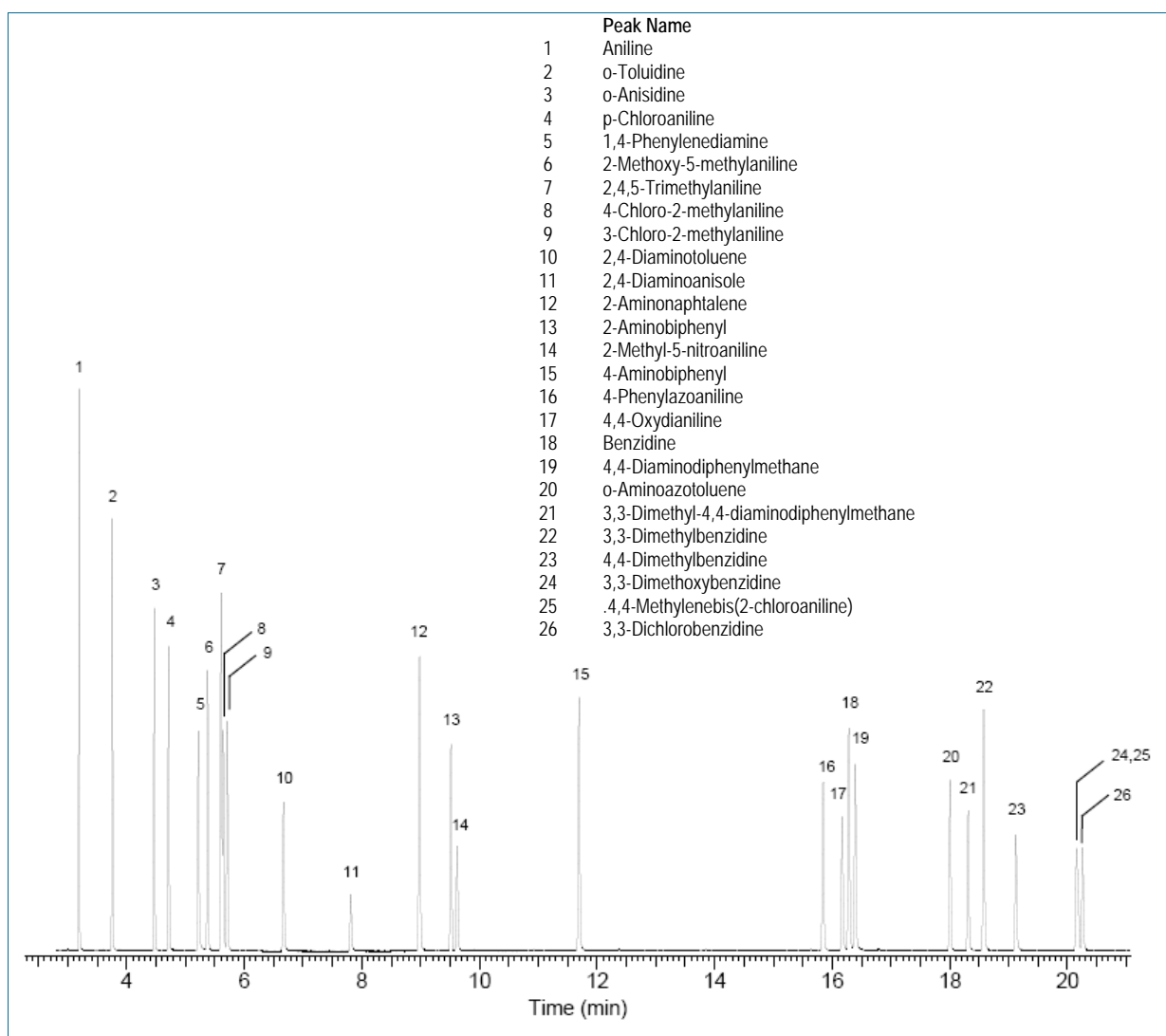

TKG 1280