


Column: **MetAmine-VOL**, P/N TR-635063
Dimensions: 60 m x 0.32 mm
Injection: 100 µL Head Space (2t,75°) split 1:15, 180°C
Sample: mix of amines in water
Carrier gas: He, 14 psi (96.5 kPa).
Program temperature: 40°C (10 min) @ 10°C/min. to 200°C
Detector: FID, 225°C

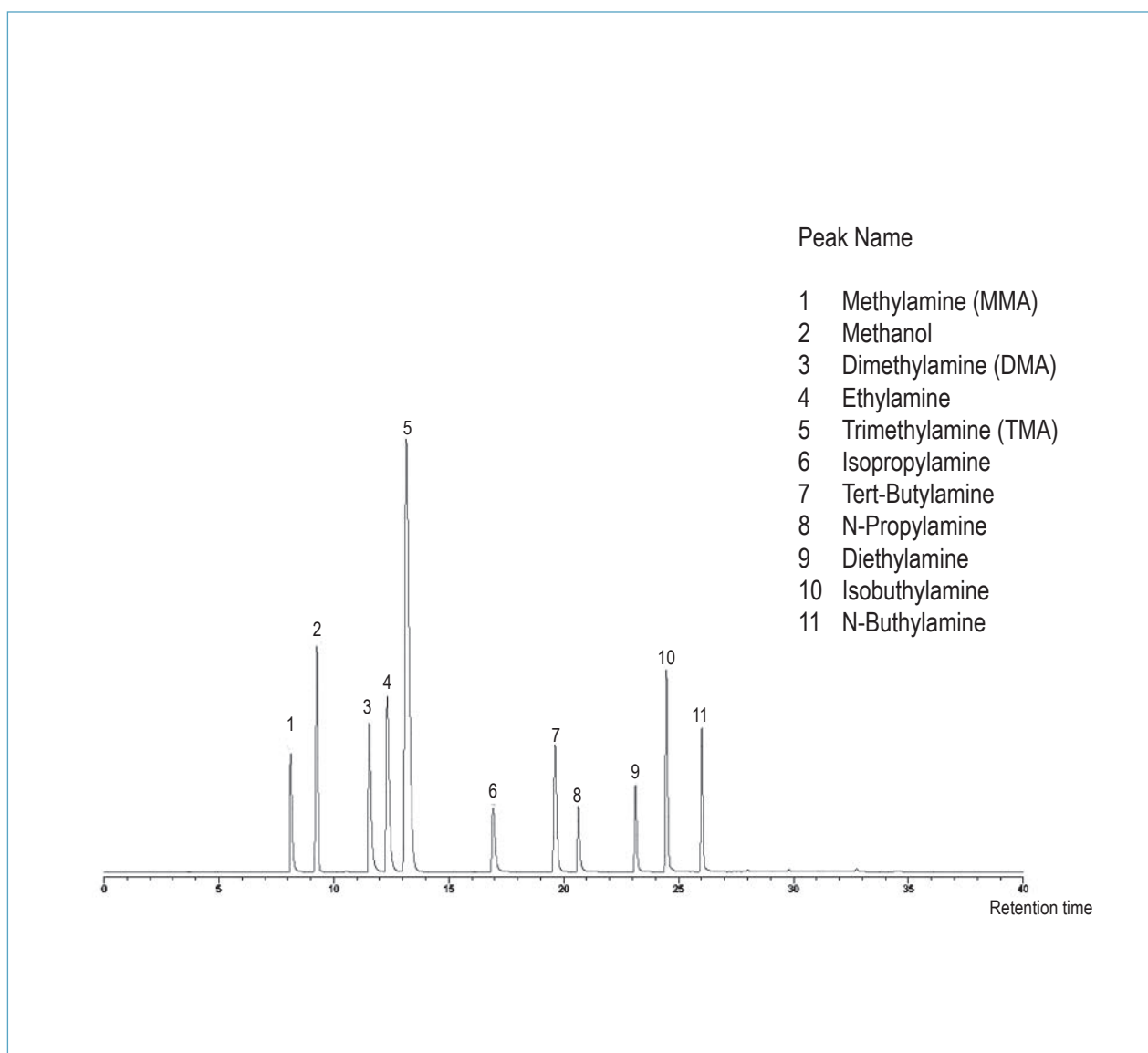

TKG 1256

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es