


**USP SOLVENTS <USP> COLUMN
TRB- G27+GUARDCOLUMN 5M**

Column: **TRB-G27**, P/N-175035
 Dimensions: 30m x 0.53mm x 5.0µm
 Oven temp.: 35°C(5') to 175°C@ 8°C/min. to 260°C (16')@35°C/min.
 Carrier gas: He, 4.5 psi (31 KPa), 35 cms. to 35°C
 Injector temp: 70°C
 FID temp: 260°C
 Injection: Direct injection of 1 µl (Uniliner), standard disolution in distilledwater (1:10)

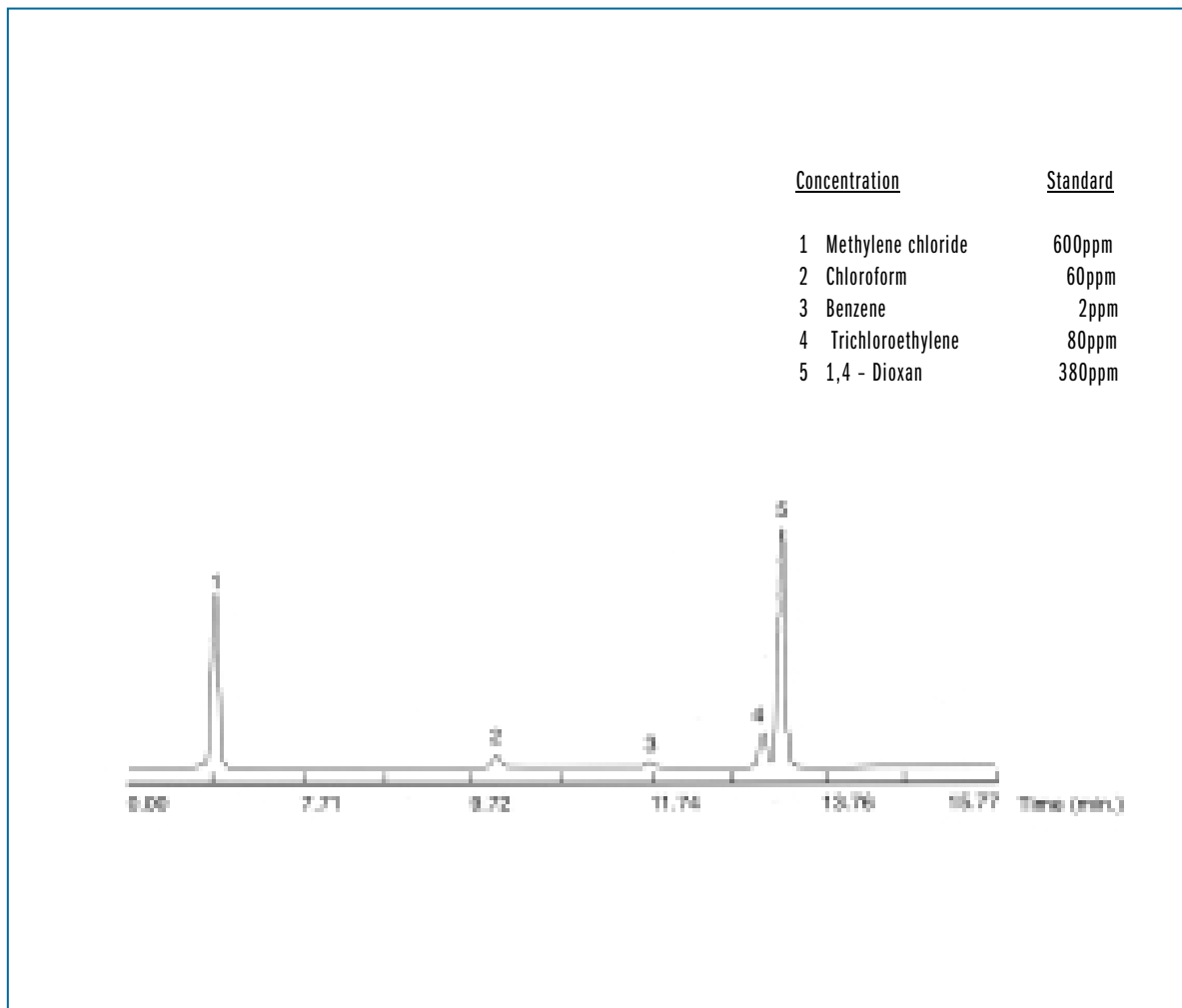

TKG 1076