

ALDEHYDES

Column: Meta. X5, P/N TR-820232


Dimensions: 30m x 0.25mm x 0.25µm

Injection: 1 µL reaction Epoxydecenal and Butylamine, splitless (1.5min), 280°C

Carrier gas: He, constant pressure 12 psi (82.7 KPa).

Oven temperature: 50°C(1min) @ 5°C/min to 240°C @ 10°C/min to 300°C

Detector: FID, 260°C

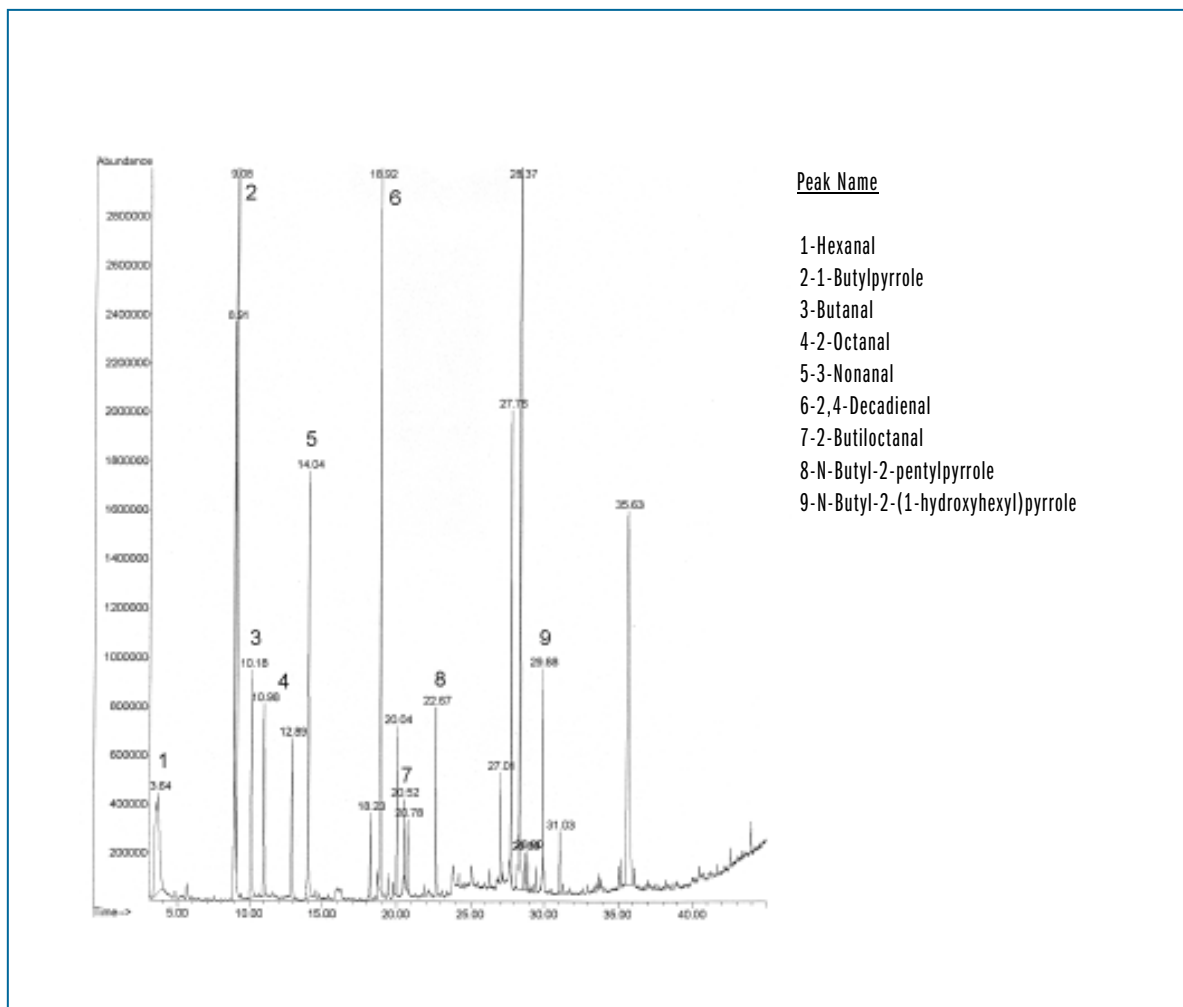

TKG 1051

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es