

## POLYAROMATIC HYDROCARBONS

Column: Meta .X5, P/N TR-821326


Dimensions: 25m x 0.15mm x 0.15µm

Injection: 0.1µL Polyaromatic hydrocarbons, (200ng/comp), splitless 30s, 300°C

Carrier gas: H2, constant pressure 35 psi (241.1 Kpa)

Oven temperature: 40°C(2min) @ 20°C/min to 200°C @ 4°C/min to 310°C(5min)

Detector: FID, 340°C

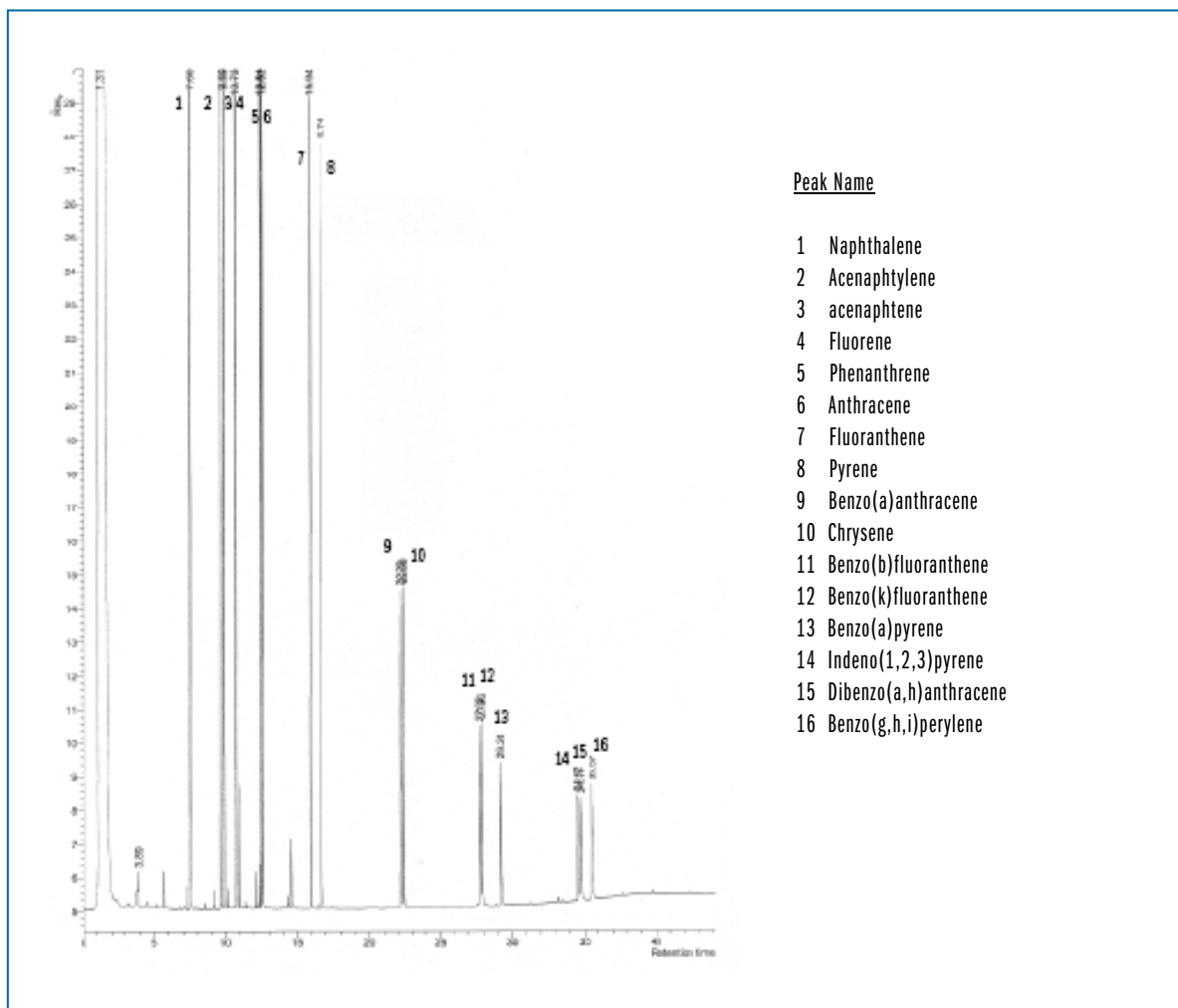

TKG 1061