


PURIFIER SLUDGE ANALYSIS

Column: **Meta X5** P/N TR-820232
Dimensions: 30m x 0.25mm x 0.25µm
Injection: 2.0 µl standard (split 1:50), 280°C
Carrier gas: He, 9 psi (62 kPa)
Oven temperature: 120°C (1min.) to 300°C (21min.) @ 10°C/min.
Detector: MS, full scan 50-550 amu, 280°C
Sample: Urban purifier sludge (250ppm di (2-ethylhexyl) phthalate, DEHP)

Chromatogram from B. Bagó, J. Díaz. Chromatography Dep. IQS.

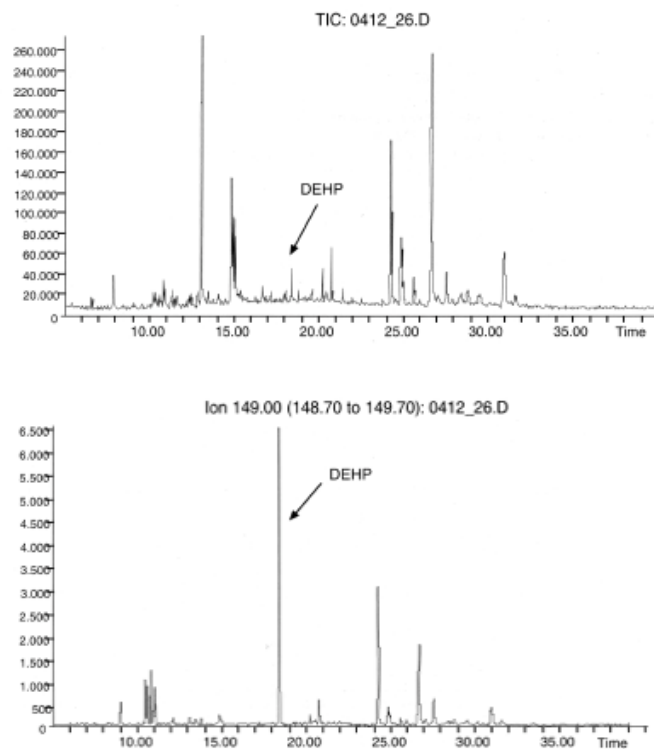

TKG 1080-1

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es