


SEPARATION OF PESTICIDES

Column: Meta .X5, P/N TR-820232
Dimensions: 30m x 0.25mm x 0.25µm
Injection: 1 µL standard (440 mg/L), split 1:20, 200°C
Carrier gas: He, 1mL/min
Oven program: 70°C(1min) @ 20°C/min to 150°C @ 3°C/min to 200°C
Detector: ECD, 250°C

Chromatogram provided by AINIA

| <u>Peak Name</u> | <u>RT (min)</u> |
|----------------------------------|-----------------|
| 1 Nonadecane (Internal Standard) | 7.91 |
| 2 Pendimethalyn | 10.09 |

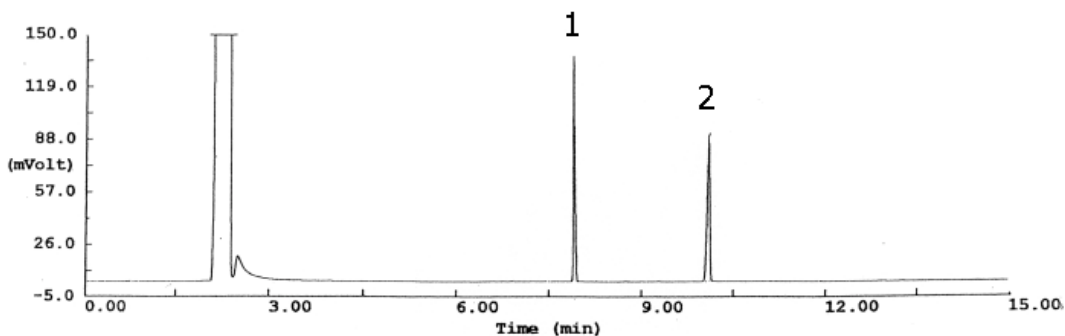

TKG 1086