


## PHENOLS-ANSIOLES IN WINE (FAST CHROMATOGRAPHY)

Column: META.X5 P/N TR-820441  
Dimensions: 10m x 0.10mm x 0.40µm  
Injection: 330°C, 0,5µL (St 10 ppb), split (200:1)  
Carrier gas: H<sub>2</sub>, 40 psi (275.6 kPa)  
Oven temperature: 100°C to 250°C @ 20°C/min.  
Detector: ECD, 330°C

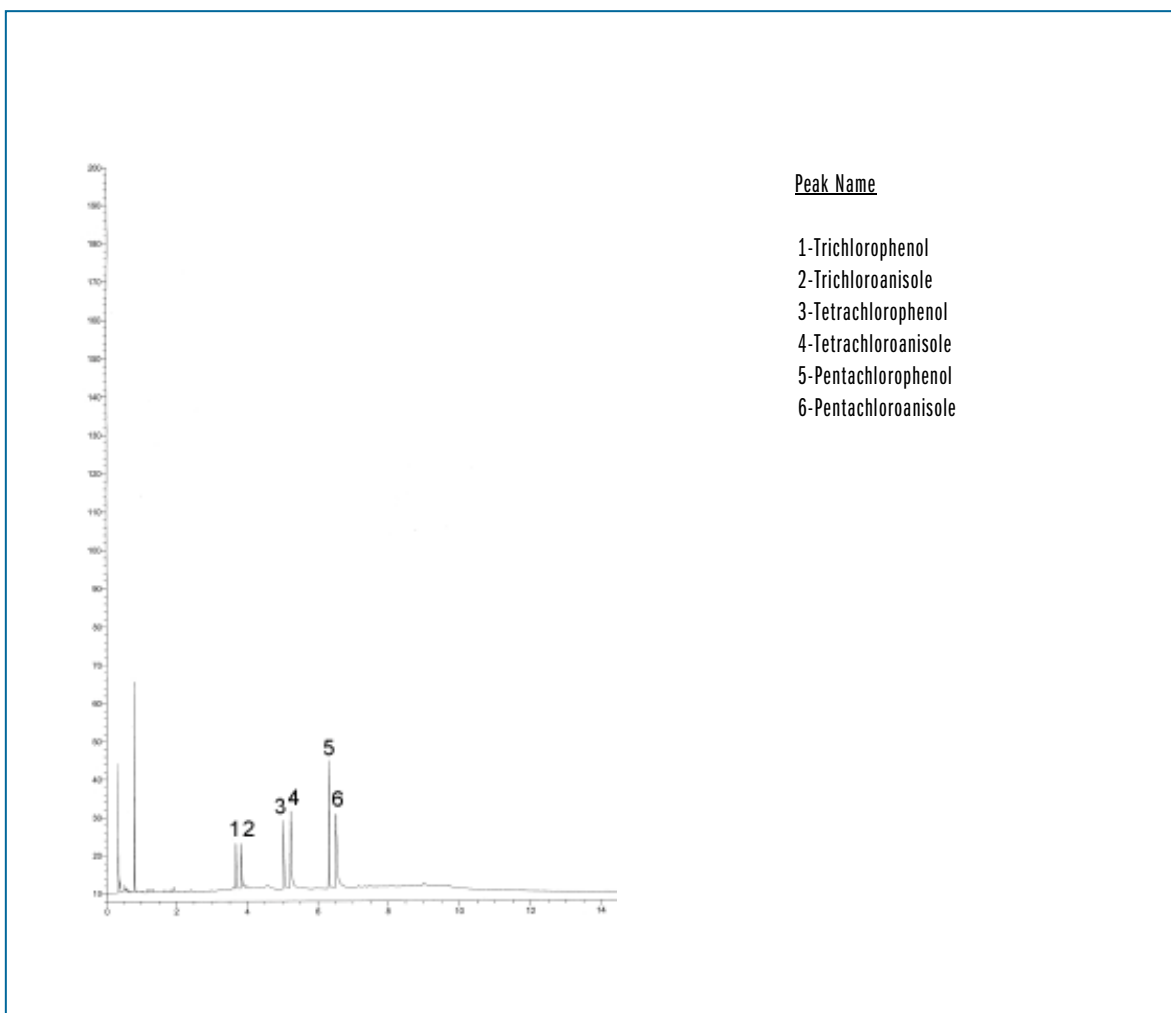

TKG 1193

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es