

SEPARATION OF PHTHALATES

Column: **Meta.X5**, 30m x 0.25mm x 0.25 μ m (P/N: TR-820232)


Carrier gas: Helium, 12psi, constant pressure mode

Injection: 250°C, split ratio 35:1

Oven temperature: 40°C (1min) to 325°C @ 12°C/min

Detector: FID, 330°C

Sample: 1 μ l Phthalate Ester Mix 1000 μ g/ml each compound

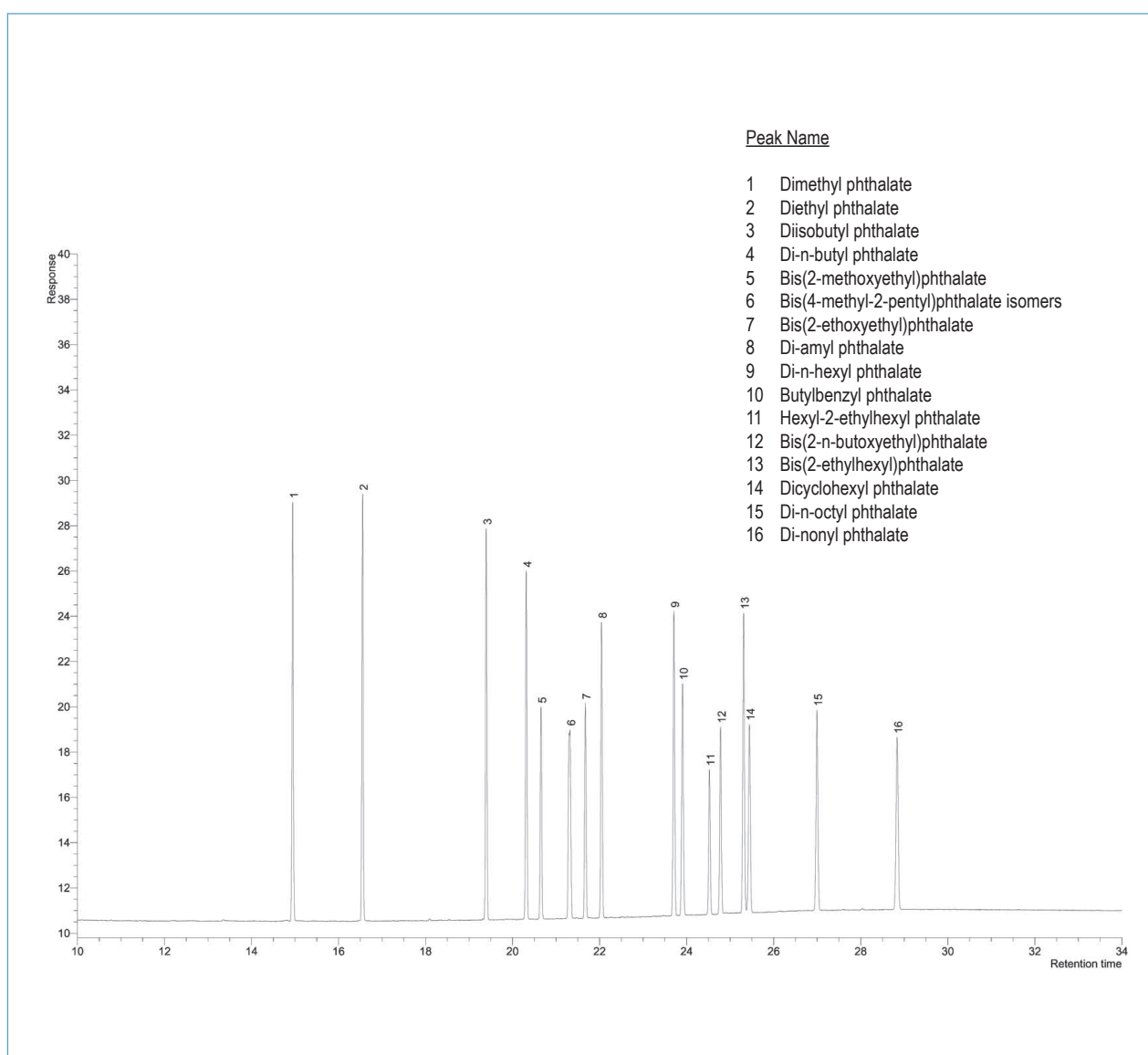

TKG 1243

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es