


SOLVENTS

Column: **TRB-624**, P/N TR-603035
Dimensions: 30m x 0.53mm x 3.0µm
Injection: wet needle, split 1:5, 260°C
Carrier gas: He, constant pressure 3 psi(20.7 KPa).
Oven temperature: 40°C (isothermal)
Detector: FID, 260°C

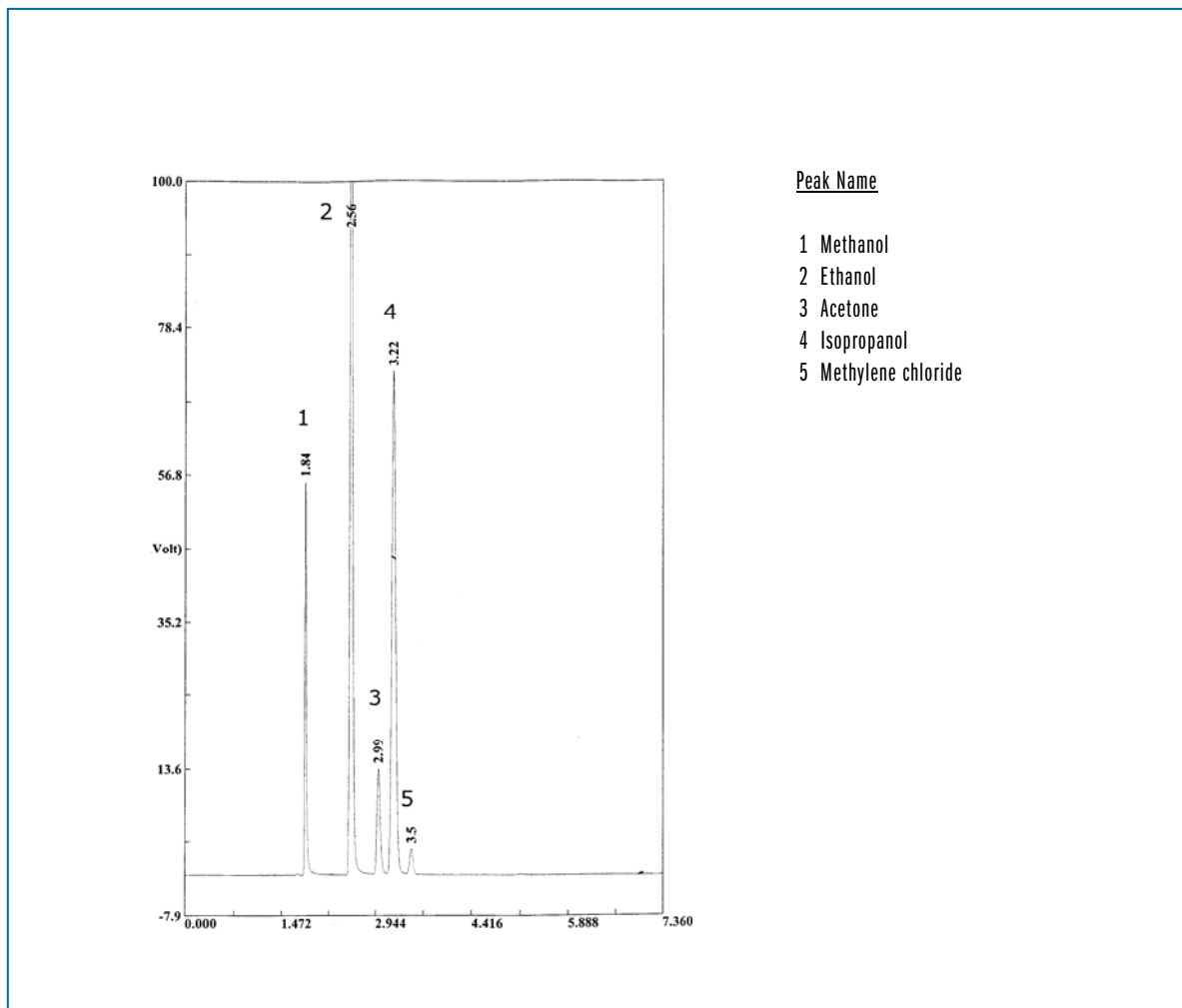

TKG 1004

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es